Before operating this unit, please read these instructions completely and save them for future use.

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# Panasonic®

Operating Instructions

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(Household) Battery-Powered Toothbrush

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Model No.

型號 EW-DS13

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
**English**

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








**24**

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**14**

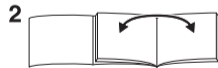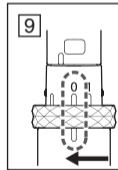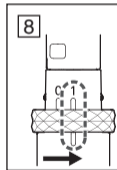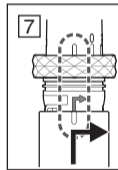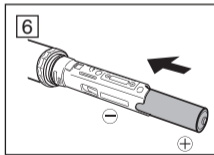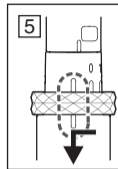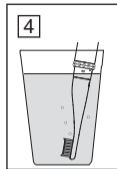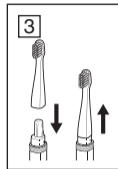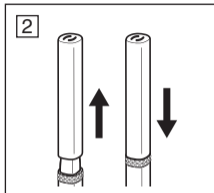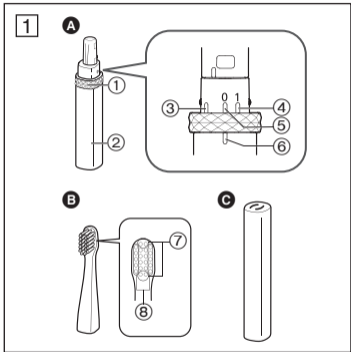Thank you for choosing a Panasonic Battery-Powered Toothbrush. Please read all instructions before use.

## Safety Precautions

### WARNING

- Keep out of the reach of toddlers and infants.
  - Failure to do so may result in accidents or accidental ingestion of the accessories or removable small parts.
- Terminate the usage when there is abnormality or

failure.

- Continued use of the toothbrush may cause fire, electric shock or other injuries.

### CAUTION

- People who have severe periodontal disease should consult a doctor before use. Failure to do so may cause damage to teeth and gums.
- Do not use the brush in combination with toothpaste that contains a lot of abrasive.
- Do not force the bristles hard against the teeth or gums, or apply the bristles to one area for prolonged periods.
- Do not apply the solid part of the brush head to teeth.

- Toddlers, infants, people who are unable to operate and people with weak oral senses should not use this appliance.
- If you have oral health problems such as teeth undergoing treatment or unstable dentures, avoid using the brush on those parts as doing so may cause damage to teeth and gums.
- Never modify, disassemble, or repair the appliance yourself as doing so may cause unexpected accidents or trouble.
- Use a mild detergent for cleaning.
- Never use thinners, benzine, or alcohol to clean the toothbrush since it may cause malfunction or cracking/discolouration of the parts.

- A guardian should provide guidance to a child on how to use the product. Make sure that the appliance is used as instructed. Failure to do so may cause damage to teeth or gums.

## Parts identification (Fig. 1)

- A** Main unit
  - ① Ring
  - ② Battery cover
  - ③ Battery cover attaching position
  - ④ Power “1” position
  - ⑤ Power “0” position
  - ⑥ Mark position
- B** Dual-Edged Brush (x1)
  - ⑦ Inverted V-shaped bristles for interdentium
  - ⑧ Extra fine bristles for between teeth and gums
- C** Cap

### Attaching and removing the cap (Fig. 2)

Attach and remove the cap straightly.

### Attaching and removing the brush head (Fig. 3)

Attach and remove the brush head straightly.

### If the brush head is hard to remove (Fig. 4)

If the brush head is hard to remove due to hardened toothpaste, soak the entire brush head in cold or lukewarm water for approximately 5 minutes and then remove.

## Replacing the battery

When the battery is low, there is less vibration and a weaker brushed feeling. Please replace the alkaline battery with a new one.

1. Remove the cap.
2. Remove the battery cover by rotating it. (Fig. 5)
  - Rotate the battery cover in the direction of the arrow, and align the marks.
  - Pull to remove the battery cover.
3. Insert an alkaline AAA (LR03) alkaline battery. (Fig. 6)

- It is possible to use the appliance for approximately 90 days by using a Panasonic alkaline AAA (LR03) alkaline battery (using the appliance once a day for 2 minutes).
4. Attach the battery cover by rotating it. (Fig. 7)
    - Align the marks and insert the main body into the battery cover.
    - Rotate the battery cover in the direction of the arrow.

**Notes**

- Batteries are not supplied. Please purchase separately.
- Before replacing the alkaline battery, ensure the unit is thoroughly dry to prevent water drops entering the main unit.
- Always use alkaline batteries. Use of non-alkaline batteries may result in a drastically reduced operating time.
- Always follow the cautions printed on the battery.
- Use only AAA (LR03) alkaline battery.
- Check the polarities (+ and -) of the battery and install it correctly.
- Remove the battery when not using the appliance for an extended period of time.
- Please dispose of used battery appropriately.
- Do not disassemble the battery or throw it into a fire. The supply terminals are not to be short-circuited.
- Remove the battery immediately and dispose of it safely when it is discharged.
- Non-rechargeable battery is not to be recharged.

## Using the toothbrush

1. Remove the cap.
2. Confirm that the brush head is fixed firmly.
3. Place the toothbrush in your mouth.
  - When using the toothbrush, especially with toothpaste applied, place the brush head in the mouth before operating to prevent splattering.
4. Rotate the battery cover to turn on the toothbrush. (Fig. 8)
5. Brush your teeth by moving the main unit by hand.
  - Brush by moving the toothbrush in the same way as when you brush by hand. Just applying the toothbrush to teeth cannot fully remove plaque, so it is recommended that you move the toothbrush.
  - Do not apply the brush head to teeth or gums too strongly.
6. When finished brushing, rotate the battery cover to turn off the toothbrush. (Fig. 9)
  - In order to prevent splattering of saliva, etc., turn off the toothbrush while it is still in the mouth.

## Tips on usage

### Boundary between the tooth surfaces or teeth and gums

- Press the bristles at a 45° angle.
- Place the brush so that the inverted V-shaped bristles can fit to the hard-to-reach spaces between teeth.
- Place the brush so that the sides of bristles can fit in between the teeth and gums.

## Tips on usage


### At the back of the front teeth or misaligned teeth

- Press the brush vertically.
- Place the brush according to the angle of the teeth.

### Bite (back teeth)

- Apply the brush to the teeth at an angle of 90° from the top.
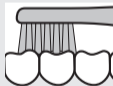
### Notes


- Do not use for any purpose other than dental care.
- Do not drop or handle roughly.
- Please be careful not to get hair products, hand cream, etc. on the main unit because it may ruin the appearance of the toothbrush.

## After use

### Washing the brush head

Wash the brush head under running water.


- Do not wash with water hotter than 80 °C.

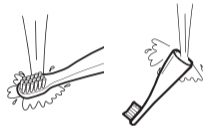- Do not pull the bristles of the brush.  
(Doing so may cause bristles to come out.)
- Shake well to remove water after washing.

**Washing the main unit**

1. Make sure that the battery cover is closed.
2. Wash the main unit under running water.

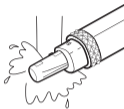- Do not wash with water hotter than 40 °C.
- Do not immerse the main unit under water. (Otherwise water may get inside.)

**Washing the cap**

Wash the Cap under running water.

- Do not wash with water hotter than 40 °C.

**Note**

- Shake the cap and brush head well to drain off water, and dry them using a dry towel, etc. After this, put the cap on the main unit. Leaving them wet may dirty or damage the pouch or bag.

## Changing the brush head (Sold separately)

Replacing the brush head every 3 months is recommended.

If the bristles are bent to one side or worn, the toothbrush is not effective for brushing.

| Part number          | Replacement Brush for EW-DS13   | Quantity |
|----------------------|---|----------|
| WEW0958-W<br>(White) |  | 2 pcs    |

- Only Replacement Brush for EW-DS13 can be used.

## Specifications

|              |  |
|--------------|--|
| Model No.    | EW-DS13  |
| Power source | AAA (LR03) alkaline battery x 1<br>(sold separately) |

This product is intended for household use only.

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© EN, TW, KR
Panasonic Corporation