
Operating Instructions
(Household use) Special Multi Treatment Device VITALIFT RF
Model No. EH-SR86

English..... 3

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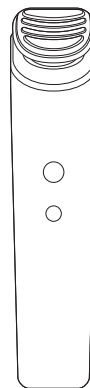Operating Instructions
(Household use) Special Multi Treatment Device VITALIFT RF
Model No. EH-SR86

Contents

Safety precautions.....	4	FAQ	29
Important information	10	Troubleshooting	32
Parts identification	11	Battery life	35
Before use	12	Removing the built-in rechargeable battery	35
Preparation before use	15	Specifications	36
How to use	17	Optional parts (Sold separately)	36
After use	28		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

The USB power adaptor is sold separately.

Safety precautions

Make sure to follow these instructions.

To prevent accidents, injuries or damage to property, please follow the instructions below.

- **The following chart indicates the degree of damage caused by improper operation.**

 **DANGER** Indicates a potential hazard that will result in serious injury or death.

 **WARNING** Indicates a potential hazard that could result in serious injury or death.

 **CAUTION** Indicates a potential hazard that could result in minor injury or property damage.

- **The symbols are classified and explained as follows.**

This symbol indicates prohibition.

This symbol indicates requirement that must be followed.

DANGER

► Preventing accidents

Do not use this appliance near electronic medical devices or home electrical therapy devices such as the following:

- **Pacemakers and other electronic medical devices implanted in the body**
- **Heart-lung machines and other life-supporting electronic medical devices**
- **Wearable electronic medical devices, such as electrocardiographs and intravenous drip equipment**

- Accidents or health problems can be caused by malfunction of the electronic medical device or electrical home therapy device.

⚠ WARNING

► To prevent skin or physical problems

The following persons should not use this appliance:

- Individuals with heart disease
- Individuals with facial nerve disorders
- Individuals with hemorrhagic disease, hemophilia or other blood disorder
- Babies, toddlers or other child in the growth stage, or other individual who cannot indicate intent
- Individuals who suffer a medical condition in the head due to a traffic accident or other traumatic event

Do not use this appliance in the RF×EMS LIFT mode near the heart.

- Doing so may cause accidents or illness.

The following persons should not use this appliance:

- Persons with allergies
- Persons with skin diseases or atopic dermatitis

- Doing so may cause skin problems.

Individuals without sense of heat should not use this appliance.

- Doing so may cause burns.


Women who are pregnant or have gone through childbirth in the past one month should not use this appliance.

- Doing so may cause skin problems due to hormonal imbalance.

Do not use this appliance in the following body area or following cases.

- When you do not feel well
 - After taking a steroid hormone, an immunosuppressant or any other medication that induces sleepiness
 - When you are very tired
 - Areas with skin damage, eczema or swelling
 - When there are abnormalities on the skin such as dermatitis or extreme sunburns
 - Body parts where epilation was performed on the same day
 - Areas of the body that have received a surgical treatment (reconstructive or cosmetic surgery)
 - Throat or mucous membranes
- Doing so may cause skin or physical problems.

<Prohibited areas>

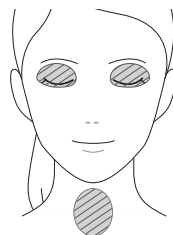Do not use the appliance on areas marked by .

Do not use on eyelids.

- Doing so may cause eye problems.

⚠ WARNING

Do not keep the appliance on the same spot of skin. Keep moving the appliance so as to continuously shift to different spots.

- Failure to do so may cause burns.

The following persons should consult own physician before using this appliance:

- Individuals with serious disease
 - Individuals with acute disease
 - Individuals with malignant tumors
 - Individuals with sensory disorders due to advanced peripheral circulatory impairments, for example those caused by diabetes
 - Individuals with infectious disease
 - Individuals with febrile disease
- Failure to do so may cause accidents or illness.

The following persons should consult own physician before using this appliance:

- Persons who have received an aesthetic or medical treatment within approximately 2 weeks that involves some downtime (the length of time that the skin needs to heal and recover) or significant reddening of the skin
 - Persons with sensitive skin
- Failure to do so may cause skin problems.

If you feel pain or abnormality on your skin or other part of your body, immediately stop using this product.

- Failure to do so may cause accidents or skin/physical problems.

**► USB power adaptor and USB cable:**

Do not immerse them in water or wash with water.

Do not place them over or near water filled sink or bathtub.

Never use the appliance if they are damaged or if the power plug fits loosely in a household outlet.

Do not damage them.

- Do not deface, modify, forcefully bend or pull them. Also, do not place anything heavy on, twist or pinch them.

- Doing so may cause electric shock or fire due to a short circuit.

Do not use them to charge other devices.

[Note that USB power adaptor (ES-RU1-02) can be used with compatible Panasonic products.]

- Doing so may cause burns or fire due to a short circuit.

Always plug the USB cable to a power supply with a rated voltage of 5.0 V and a rated output current of 1.0 A which meets requirements of the Safety Extra Low Voltage (SELV) system.

- Failure to do so may cause burns or fire due to a short circuit.

When cleaning the appliance, disconnect the USB power adaptor from the household outlet.

- Failure to do so may cause electric shock or injury.

⚠ WARNING

Regularly clean them to prevent dust from accumulating.

-  - Failure to do so may cause fire due to insulation failure caused by humidity.
Remove the USB power adaptor and wipe it with a dry cloth.

Use the specified USB power adaptor and USB cable.

-  - Failure to do so may cause burns, electric shock or fire due to a short circuit.

► Power supply

Do not connect or disconnect the USB power adaptor from the household outlet with wet hands.

-  - Doing so may cause electric shock or injury.

Do not use in a way that would exceed the rating of the household outlet or the wiring.

-  - Exceeding the rating by connecting too many plugs to one household outlet may cause fire due to overheating.

Fully insert the power plug of the USB power adaptor securely.

-  - Failure to do so may cause fire or electric shock.

► This product

This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, leave, or store in a high temperature environment.

- Doing so may cause overheating, ignition, or explosion.

Do not modify or repair.

- Doing so may cause fire, electric shock, or injury.
Contact an authorized service center for repair (replacement of battery etc).

Never disassemble except when disposing of the product.

- Doing so may cause fire, electric shock, or injury.

► In case of an abnormality or malfunction

Immediately stop using the device and disconnect the USB power adaptor.

- Failure to do so may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

- The main unit, USB power adaptor or USB cable is deformed or abnormally hot.
- The main unit, USB power adaptor or USB cable smells of burning.
- There is an abnormal sound during use or charging of the main unit.
- Immediately request inspection or repair at an authorized service center.

⚠ WARNING**► Take care of the following precautions as well**

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or

instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

- Failure to do so may cause an accident or injury.

⚠ CAUTION**► Power supply**

Do not charge when the appliance socket is wet.

- Doing so may cause burns or fire.

Be careful that foreign objects such as dust or metallic pins do not get on the appliance's power plug, USB plug, USB port and appliance socket.

- Doing so may cause electric shock or fire due to a short circuit.

Disconnect the USB power adaptor from the household outlet when not charging.

- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

► Note the following precautions

Do not apply the appliance's head to a metallic surface.

- Doing so may cause burns.

Before using, check that the head is free from cracks, chips, deformities, and other damage.

- Failure to do so may cause injury to the skin.

Use only for care of the face, chest and neck (excluding the arm used for a patch test).

- Failure to do so may cause skin damage, inflammation, etc.

If you feel discomfort in your teeth while using the appliance, select a lower treatment level or avoid using the appliance on the discomfort-causing spot of the skin.

- Use of the appliance may cause temporary tooth pain.

Comply with the recommended use frequency and the manipulation method specified in this manual.

- Failure to do so can damage the skin or can cause or aggravate other problems such as infection, inflammation, bleeding or abnormal pigmentation of the skin.

Handling of the removed battery when disposing

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- **Never peel off the battery cover film.**
 - Doing so may cause overheating, ignition, or explosion.

WARNING

After removing the rechargeable battery, do not keep it within reach of children and infants.

- The battery will harm the body if accidentally swallowed.
If this happens, consult a physician immediately.

If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.

- Failure to do so may cause eye injury.
Consult a physician immediately.

CAUTION

If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.

- Failure to do so may cause inflammation.

Important information


- Charge the appliance correctly according to these Operating Instructions. (See page 16 "Charging".)
- This appliance contains battery that is only replaceable by skilled persons. Please contact an authorized service center for the details of repair.

Note the following precautions in order to prevent malfunction and damage

- Do not drop or hit against other objects.
- Store the appliance in a place with low humidity after use.

This symbol on the product means that not generated higher voltage than the safety extra low voltage (SELV). In which protection against electric shock relies on power supply.

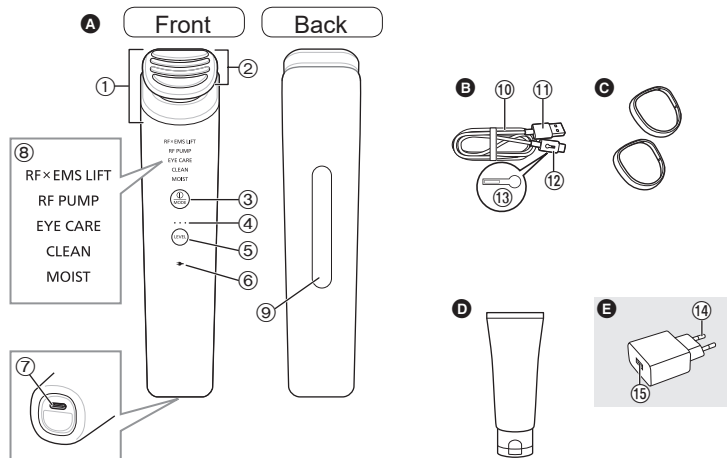A Main body

- ① Type D lift head (the "head")
- ② Electrodes
- ③ POWER/MODE switch
- ④ Level indicator [white]
- ⑤ Level selector switch
- ⑥ Charging indication/notification lamp (🔌)
- ⑦ Appliance socket (Type-C)
- ⑧ Mode indicator lamp [white]

⑨ Ion inductive panel

- Senses that the Ion inductive panel and the head are in contact with the skin and outputs EMS, RF and ion.

B USB cable

- ⑩ Cord
- ⑪ USB plug (Type-A)
- ⑫ USB plug (Type-C)
- With overheat protection function
- ⑬ Overheat protection function mark

Accessories

C Cotton ring (x2)

D Included gel

Conductive Gel (EH-4R03)

Optional parts

E USB power adaptor (ES-RU1-02) (Sold separately)

- ⑭ Power plug
- ⑮ USB port (Type-A)

Before use

► Check your skin

- Before use, a skin allergy test is required. (See page 15.)
 - This quick test checks for any potential skin issues before you begin using the appliance normally.
- Please use this appliance with bare skin (after removing make-up and washing your face).
- Excessive use of the appliance can be detrimental to the skin. Monitor your skin condition and adjust the amount of use accordingly.
 - Once any problem is found with the skin, immediately stop using the appliance and consult a physician.
 - Skin problems can be caused by various factors such as your personal physical properties or health condition, mental or emotional stress, internal organ disorders or other beauty products you are using. It is important to consult a medical specialist to properly address skin problems.

► Treatment levels

- Select the desired level according to the skin conditions and personal preferences.
 - The choice of treatment levels is different depending on which treatment mode you are using. (See page 13.)
 - For the first use, it is recommended to set the level to "1".

► Sensors

- This appliance is equipped with various sensors to assist safe use.

- Overuse avoidance sensor

If the head is kept applied on the skin for more than 8 seconds, this sensor temporarily shuts down the appliance's RF and/or EMS output. The output will resume after the appliance is moved away from the skin once and applied again.

- Temperature sensor

When the temperature gets too high, this sensor temporarily stops the RF output for safety. The output will resume when the head temperature drops down to a safe level.

► Accessories

- Before use, remove earrings or any other jewelry from the face or around the neck.

► Use of the included gel or other skin care products

- This appliance can be used with the included gel or other skin care products that you prefer.
 - Use the appliance with a product that reduces skin stress and assists smooth sliding of the head over the skin to help deliver EMS and RF to the skin most effectively.

Notes

Do not use the appliance with the following types of products:

- Oil-based products, products that contain scrubbing particles, lotion with a warming sensation

Treatment modes

► Intensive treatment

Mode	Features	Technologies employed	Treatment time
RF×EMS LIFT	To exercise facial muscles	EMS*, RF, ion and LED	Approx. 6 minutes
RF PUMP	Warming treatment with RF	RF*, ion and LED	Approx. 6 minutes
EYE CARE	Skin care around the eyes	Ion*, RF and LED	Approx. 2 minutes

► Daily treatment

Mode	Features	Technologies employed	Treatment time
CLEAN	To remove dead cell buildup from the skin surface and to cleanse skin pores	Ion*, RF and LED	Approx. 3 minutes
MOIST	For penetrative skin treatment through the stratum corneum layer	Ion*, RF and LED	Approx. 3 minutes

For each feature marked with “*”, you can select one of 3 treatment levels by using the level selector switch.

Use frequency

The allowable frequency of use differs depending on the treatment mode.
Comply with the recommended use frequency and refer to the example use patterns shown below.
Do not use the appliance in an intensive treatment mode more than once in a single day.

► Intensive treatment

Mode	Use frequency	Weekly use pattern example						
		Mon	Tue	Wed	Thu	Fri	Sat	Sun
RF×EMS LIFT	Select one or the other mode once a day • Use up to 3 times per week in 2 modes total. • Wait at least 1 day between each use.		✓					✓
RF PUMP						✓		
EYE CARE	Once a day • Use up to 3 times per week. • Wait at least 1 day between each use.	✓			✓		✓	

(✓: Use)

► Daily treatment


Mode	Use frequency	Weekly use pattern example						
		Mon	Tue	Wed	Thu	Fri	Sat	Sun
CLEAN	Once a day • Each mode can be used everyday • Both of the modes can be used on the same day in combination with an intensive treatment mode.	✓	✓	✓	✓	✓	✓	✓
MOIST		✓	✓	✓	✓	✓	✓	✓

(✓: Use)

Preparation before use

Skin (allergy) test procedure

To avoid potential skin problems that may occur when the included gel or other skin care product is used in combination with this appliance, be sure to run a patch test before full use.

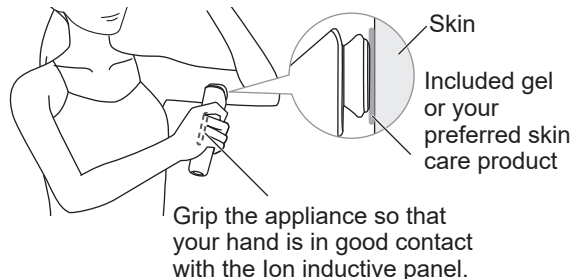- 1 Wash the area of skin to be tested (inside of the arm) with soap and wipe the skin dry with a towel.**

- 2 Apply the included gel or other skin care product you want to test.**

- Apply to make a thin film about the size of the appliance's head.
- If you want to test the included gel and another skin care product at the same time, apply each on the skin with some distance apart.

Included gel

Your preferred skin care product

- 3 Turn the appliance power on and use it on the area of the skin where the gel or skin care product is applied. (See page 21.)**

1. Turn the appliance power on.
 - If testing with the included gel: Select the RF×EMS LIFT mode and level “1”.
 - If testing with other skin care products: Select the RF PUMP mode and level “1”.
2. Apply the head to the skin and move in a slow circular motion (approx. 7 seconds).
 - If you are testing the included gel and another skin care product at the same time, wipe the head clean after testing with the first product before testing with the other one.
 - After finishing the test, let the skin naturally dry.

- 4 Observe the skin for 48 hours after the test.**

- Do not rub the area of skin being tested.

► If you notice a skin problem within 48 hours:

Do not use the appliance in combination with the included gel or the skin care product that you have used in the test.

Charging


Wipe the appliance socket dry if you see beads of water or skin care product residue in or around the appliance socket.

Use the specified USB power adaptor and USB cable.

1 Connect the USB cable to the USB power adaptor and the appliance's main body.

2 Insert the USB power adaptor fully into the household outlet.

- Check that the charging indication/notification lamp (🔌) is glowing.
- Charging is completed after approx. 3 hours.

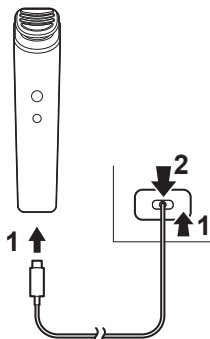3 When the charging indication/notification lamp (🔌) turns off, disconnect the USB power adaptor from the household outlet and then disconnect the USB plug from the appliance.

► To check if charging is completed

Disconnect the USB plug once and then reconnect while charging.

→ If the appliance is already fully charged, the charging indication/notification lamp (🔌) will turn on once and then turn off after approximately 5 seconds.

► When to recharge

If the charging indication/notification lamp (🔌) starts to flash (once per second), this indicates that the battery charge level is getting low. Recharge the appliance.

- Once fully charged, the appliance will be good for approximately 4 days of use. (See page 14.)
Note that the appliance may last for fewer days depending on the operating environment or usage.

Notes


- When charging the appliance for the first time or when it has not been in use for more than 6 months, the charging time may change or the charging indication/notification lamp (🔌) may not glow for a few minutes. It will eventually glow if kept connected.
- The required charge time may vary depending on the battery condition.
- Recommended ambient temperature for charging is 15 °C – 35 °C. It may take longer to charge or the battery may fail to charge properly in extremely high or low temperatures.
- If you have noise problems on your radio or audio player while using or recharging the appliance, use or charge the appliance elsewhere.

How to use


Use it appropriately depending on your skin condition and the type of skin problem you want to address. Refer to P.14 for the recommended use frequency.

Basic steps for use

► If starting from the CLEAN mode

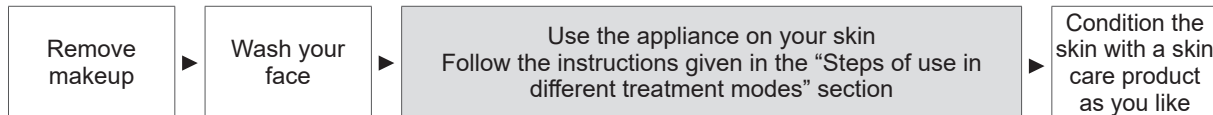► If starting from a treatment mode other than the CLEAN mode

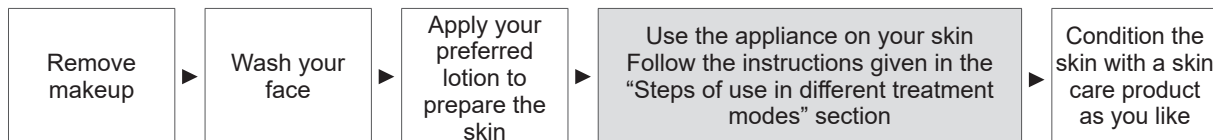Combination of multiple treatment modes

CLEAN mode and MOIST mode can be used on the same day in combination with an intensive treatment mode. For such combined use, use in the order of CLEAN mode → MOIST mode → intensive treatment mode.

Steps of use in different treatment modes

► Intensive treatment

RF×EMS LIFT mode

RF PUMP mode

EYE CARE mode

► Daily treatment

CLEAN mode


MOIST mode

*1 Recommended amount of use of the included gel (typical): 1 to 2 grams (= approximately one 1.5 cm-diameter dollop). Use of the included gel is most highly recommended as it helps deliver EMS to the skin most effectively.

<Typical recommended amount of use>

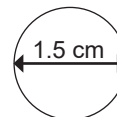
*2 The beep lets you know that half the treatment time has passed.

Preparing a cotton pad and your preferred skin care products (for CLEAN mode only)

Before using the CLEAN mode, perform the steps below:

1 Attach an ordinary makeup cotton pad to the head.

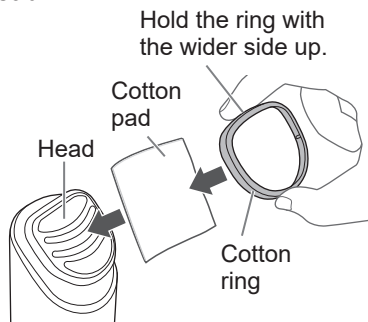Check that the cotton ring is facing the right way.

The inward-curving side of the ring should be facing outward

Press the cotton ring onto the head to engage.

Do not push in too deep.

2 Soak the cotton pad with your preferred skin lotion (apply a generous amount so that the liquid soaks through to the other side of the pad).

Using the appliance

Remove makeup and wash your face before using the appliance.


To start using the appliance from a mode other than the CLEAN mode, condition the skin with your preferred skin lotion before applying the included gel or your preferred skin milk or nutrition product.

1 Long-press the button (approx. 1 second) to turn the appliance power on.

- The [RF×EMS LIFT] lamp will come on.

2 Select the desired mode by pressing the button.


- Each time the button is pressed, the mode switches to the next one in the sequence below:

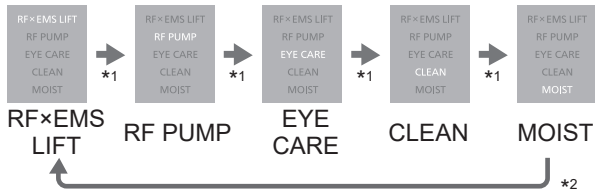[*1 Single beep/*2 Two beeps]

3 Select the desired level by pressing the button.

- Each time the button is pressed, the mode switches to the next one in the sequence below:
- Start with level 1.

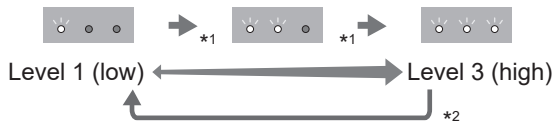[*1 Single beep/*2 Two beeps]

4 Use the appliance by moving the head over the skin

- Grip the appliance so that your hand is in good contact with the Ion inductive panel.
- Use the appliance by gently pressing the head onto the skin. Keep moving the head so as not to remain on the same spot of skin.
- Do not press the head too hard against the skin.

Ion inductive panel

Tips

When doing the cheeks and the jawline, apply the flat part of the head to the skin and move the head in an upward pulling motion.

Notes

- You cannot operate when the appliance is connected to the household outlet. Disconnect the USB power adaptor from the outlet to use the appliance.
- Use the appliance in an ambient temperature approximately 15 °C or higher. The head may not warm up when the operating temperature is not within the recommended range.

- If your hand is away from the Ion inductive panel or the head is away from the skin for more than approximately 15 seconds, a “bleep-bleep, bleep-bleep” tone will alert you. If the condition continues for approximately 30 seconds, the appliance power will automatically turn off with a “bleep-bleep, bleep-bleep and a (long) beep” tone.
- When the treatment time for the selected mode is expired, the appliance will automatically turn off with a long “beep” tone.
If you want to turn the appliance off before the treatment time expires, long-press the  button.
- After using the CLEAN mode, remove the cotton pad from the head.
- Do not reuse the removed cotton pad. (Doing so may cause skin problems.)
- If you want to switch to a different mode halfway through the treatment, turn the appliance power off once before selecting the next mode.
If you switch to a different mode without turning the appliance power off, the appliance will turn off automatically with a long “beep” tone approximately 6.5 minutes after the first power-on. This prevents the excessive use.

Tips to use EMS effectively

- “EMS” stands for “electrical muscle stimulation”, which refers to the use of electrical impulses to cause human muscles to contract. EMS stimulation differs depending on the area of the skin being treated or the type and the amount of the skin care product you are using with the appliance.

- If you are not satisfied with the perceived strength of the EMS, try the following 5 steps:

1. Condition your skin with a lotion

Wash your face and then condition the skin with the lotion before using EMS.

2. Wet the skin well

Apply a generous amount of the included gel or your preferred skin care product.

If the skin is not sufficiently wetted or starts to dry, EMS stimulation may not feel strong enough. Add the gel or skin care product as required to keep the skin wet. Use of oil-based products is not recommended as they tend to weaken the perceived strength of EMS.


Included gel

Your preferred skin care product

3. Select RF×EMS LIFT mode


This appliance outputs EMS in the RF×EMS LIFT mode only.

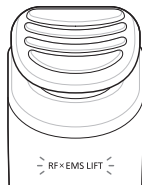
4. Make sure that the appliance is detecting the skin

1. Make sure that your hand is in contact with the Ion inductive panel.

Grip the appliance so that your hand is in good contact with the Ion inductive panel.

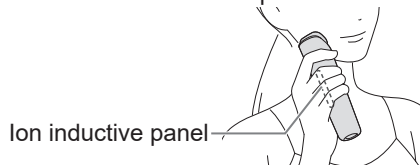2. Make sure that the electrodes touch the skin. You cannot feel EMS stimulation if the electrodes are not in close contact with the skin.

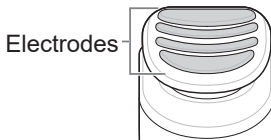5. Comply with the maximum skin contact time

If the head is kept applied on the skin for more than 8 seconds, the RF and/or EMS output will be temporarily suspended. The output will resume after the appliance is moved away from the skin once and applied again.
(Overuse avoidance sensor function)

Notes

- When the skin is dry, EMS stimulation may be harder to feel.
- If the perceived strength of EMS is not strong enough when the appliance is used with a commercial skin care product, use of the included gel is recommended.

How to apply


Treat one side of the face first before doing the other side.

- When half the treatment time for the selected mode has passed, a “beep-beep” tone will alert you. Use this as a cue to start treating the other side of the face.
- Treat each side of the face the same way by following the numbered stroke paths in the picture.

RF×EMS LIFT/RF PUMP modes

Approx. 6 minutes

Power turns off automatically

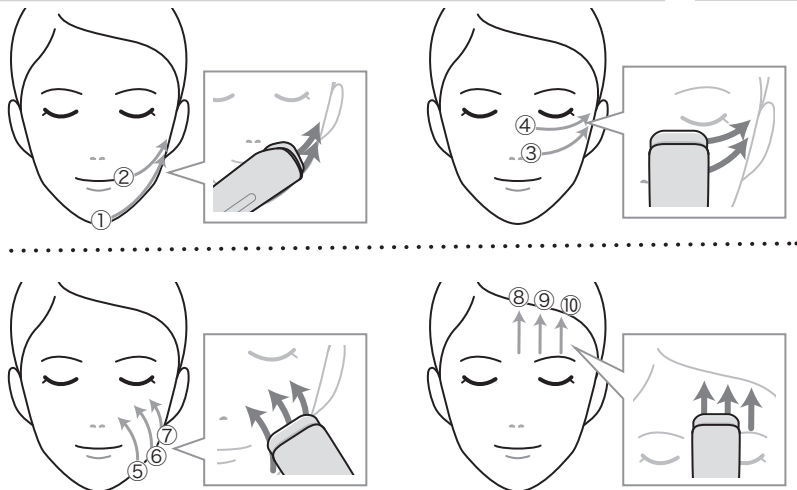
① to ⑦: Approx. 7 seconds per path, ⑧ to ⑩: Approx. 3 seconds per path

- Move the head slowly in an upward motion.
- Do not repeat doing the same path. Keep shifting to different paths in the order of ① to ⑩.

EYE CARE mode

Approx. 2 minutes

Power turns off
automatically

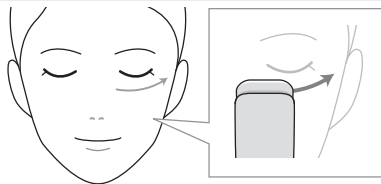
Approx. 5 seconds per numbered stroke path

- Do not use on eyelids.


CLEAN mode

Approx. 3 minutes

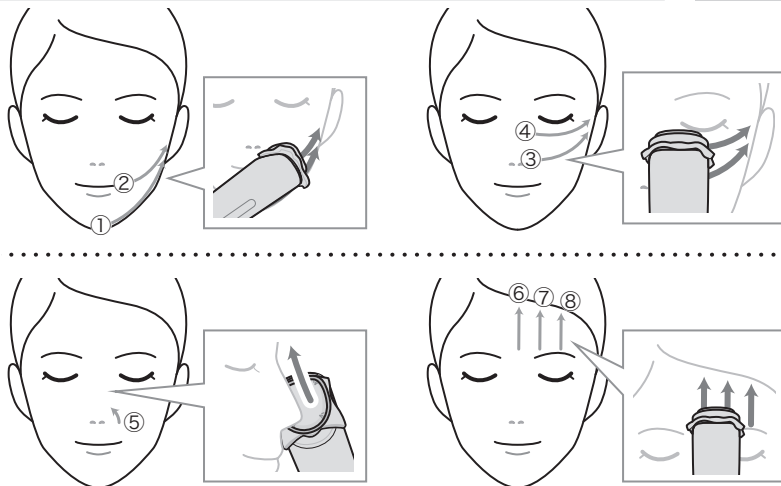Power turns off
automatically

① to ③: Approx. 7 seconds per path, ④: Approx. 5 seconds, ⑤ to ⑧: Approx. 3 seconds per path

- Do not repeat doing the same path. Keep shifting to different paths in the order of ① to ⑧.

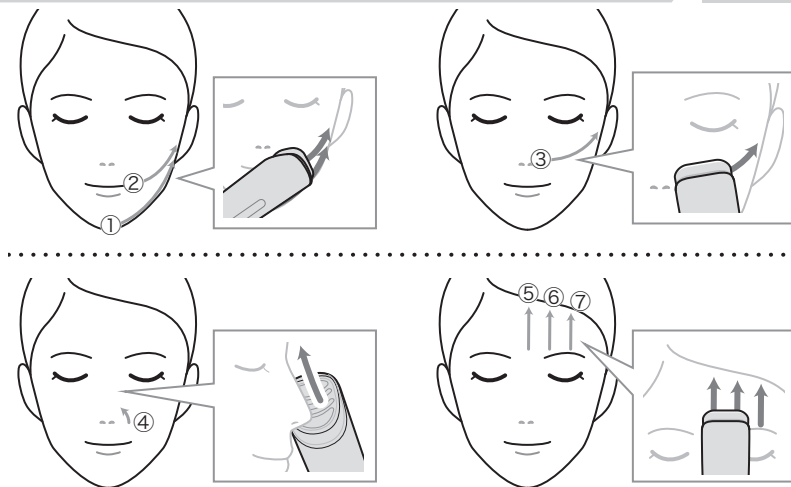① to ③: Approx. 7 seconds per path, ④ to ⑦: Approx. 3 seconds per path

- Do not repeat doing the same path. Keep shifting to different paths in the order of ① to ⑦.

Notes


- If you cannot slide the head smoothly over the skin when using the RF PUMP or EYE CARE mode, add the skin care product as required.
- If the perceived strength of EMS stimulation is not strong enough when using the RF×EMS LIFT mode, add the included gel or your preferred skin care product as required.
- If the perceived strength of EMS is not strong enough when the appliance is used with a commercial skin care product, use of the included gel is recommended.
- Make sure that your hand is in good contact with the Ion inductive panel.
- Make sure that the electrodes are in close contact with the skin. (See page 23.)
- If the head is kept applied on the skin for more than 8 seconds, the RF and/or EMS output will temporarily stop. (See page 12.)

- After using the CLEAN mode, remove the cotton pad from the head.
- Do not reuse the removed cotton pad.
- If the cotton pad starts to dry and the head does not slide smoothly while using the CLEAN mode, add the lotion to the cotton pad as required.
- If you cannot slide the head smoothly over the skin when using the MOIST mode, add the skin care product to the skin as required.

After use

Condition the skin

- Use your regular skin care product to condition the skin after using the appliance.

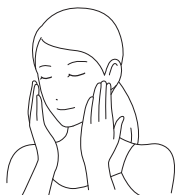Wash the appliance's main body with water

- Check that the power is off before washing the appliance.
- If the appliance is still not clean enough, dissolve some amount of ordinary bar soap in water and wipe the appliance with a cloth soaked with the soapy water. Afterward, rinse the appliance well with clean water.


1 Rinse the head and the cotton ring with water to remove residue of the gel or the skin care product.

- Do not submerge the appliance in water.



2 If there are water droplets in or around the appliance socket, remove water in the following way:

- With the appliance socket area facing down, lightly tap the appliance onto a dry cloth to drain water from the appliance socket.
- Be careful not to charge the appliance with the appliance socket wet.

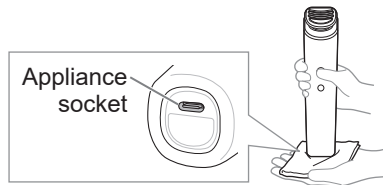3 Wipe the appliance with a dry cloth and let it thoroughly dry in a well-ventilated place.

- Do not store the appliance with the USB cable wound around the main body.

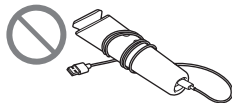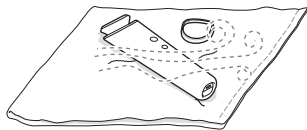Notes

- Do not use anything other than water and bar soap to clean the appliance. (For example, do not use detergents, alcohol, or nail polish removers.)
- Thoroughly remove any gel or skin care product residue from the appliance. (Failure to do so may cause malfunction or component crack/discoloration.)

FAQ

Question	Answer
What is RF?	▶ RF stands for “radio frequency”, which refers to high-frequency electromagnetic radiation. RF output causes a warming sensation in the skin.

Question	Answer
What is EMS?	▶ EMS stands for “electrical muscle stimulation”, which refers to a series of electrical impulses that causes human muscles to contract.
What is ION?	▶ “Ion” refers to a treatment technology that uses low-level electrical current that helps cleanse the skin surface and skin pores and also accelerates delivery of skin nutrients through the stratum corneum layer.
What is LED?	▶ “LED” refers to a type of phototherapy for skin treatment.
Can I use the appliance on my hands and body?	▶ Use the appliance only for beauty treatment of your face, neck and decollete (upper chest area).
Can I use the appliance as many times as I want in a single day?	▶ The allowable frequency of use differs depending on the treatment mode. (See page 14.)
Can I use the CLEAN mode to remove makeup?	▶ No, you cannot. Use it only after removing makeup and washing your face.

Question	Answer
Why do I have to apply a gel or skin care product on my skin to use the appliance?	▶ Applying a gel or skin care product helps deliver EMS or RF to the skin most effectively.
Is it OK to use other gel or skin care products when using the RF×EMS LIFT mode?	▶ Use of the included gel is best recommended to help deliver EMS stimulation to the skin most effectively, but it is also acceptable to use your preferred skin milk or nutrition product in its place. (See page 18.)
Do I have to rinse the included gel off my face after use?	▶ No, you don't. However, you can wipe off or rinse off the gel if you don't want the residue on your skin when applying other skin care products after using the appliance.
I am getting crumbly stuff on the skin when I apply foundation on my face after using the included gel.	▶ Wash the gel off the skin before applying makeup.
What kind of skin care product is recommended for use after using the appliance?	▶ Use your regular skin care item as appropriate depending on your skin condition.

Question	Answer
What if I feel irritation with the gel or skin care product when using the appliance or start having other skin problems?	▶ Do not use the appliance until your skin heals and recovers. (Consult a medical specialist.) • Such condition could be an indication of increased skin sensitivity. • Skin problems can be caused by various factors such as your personal physical properties or health condition, mental or emotional stress or internal organ disorders.
The appliance feels too hot.	▶ • Select a lower treatment level or move the appliance a little quicker over the skin. • Start using the appliance from level "1". • Do not keep the head applied to the same spot. Always keep moving the head over the skin when using the appliance.
The appliance doesn't feel warm enough (CLEAN mode)	▶ A thick cotton pad may weaken the warming sensation of RF. Try using a different cotton pad.

Question	Answer
EMS stimulation feels too weak.	<ul style="list-style-type: none"> • If the applied amount of the included gel or skin care product on the skin is not enough or the skin feels dry, add the special gel or skin care product to the skin as required. • Use of oil-based products is not recommended as they tend to weaken the perceived strength of EMS stimulation. • If the perceived strength of EMS is not strong enough when the appliance is used with other skin care products, use of the included gel is recommended. • Select the RF×EMS LIFT mode. • Make sure that your hand is in good contact with the Ion inductive panel.
	<ul style="list-style-type: none"> • Make sure that the electrodes are in close contact with the skin. • If the head is kept applied on the skin for more than 8 seconds, the RF and/or EMS output will temporarily stop. (See page 12.)

Question	Answer
EMS stimulation feels too strong.	<p>▶ Select a lower level.</p>
What if I feel tooth pain when using the appliance?	<p>Use of the appliance may sometimes irritate tooth nerves. Having a dental filling may worsen this tooth pain.</p> <p>▶ • Select a lower level. • Be careful not to press the head too strongly to the skin, and avoid using the appliance in areas that cause tooth pain.</p>
What is the service life of the battery?	<p>The service life of the product battery is approximately 3 years.</p> <p>▶ • If you start to notice that the appliance lasts for smaller number of uses even after a full recharge, this may be an indication that the battery is approaching the end of its service life. (The battery service life may vary depending on the use and storage conditions.)</p>

Question	Answer
Can I change the battery by myself?	▶ No, you cannot. (Doing so may cause malfunction.) Contact the store where you purchased the unit.
Does the battery deteriorate if I don't use the appliance for an extensive period of time?	▶ Yes, it can. If you are not using the appliance for more than 6 months, fully recharge the battery at least once every 6 months.
If I fully recharge the battery every time, does it affect the battery service life?	▶ No, it does not. As the appliance uses a lithium-ion battery, repeated full recharges do not affect the battery service life.
Is it OK to recharge the appliance with the appliance socket wet?	▶ No, it is not. Wipe the appliance to remove water or skin care product residue from the appliance socket area and allow the appliance to dry before recharging. (See pages 28 to 29.)

Question	Answer
Can I keep the appliance connected for recharge all the time when not in use?	▶ Keeping the appliance connected for recharge for more than 3 hours does not adversely affect its performance. For best safety and energy saving, however, it is recommended to disconnect the USB power adaptor after the appliance is fully charged.
My cotton pad comes off easily from the head.	▶ Some types of cotton pads do not adequately fit the appliance head due to their size, thickness or shape. Try using a different cotton pad.

Troubleshooting

Perform the following actions.

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Problem	Action
<p>The power does not turn on even after  is pressed.</p>	<p>► Check if you are using a specified USB power adaptor or USB cable. → Always use the specified power adaptor and cable.</p>
	<p>► Charge the appliance.</p>
	<p>► Check if the POWER/MODE switch is pressed long enough (approximately 1 second).</p>
<p>The charging indication/notification lamp () does not turn on or goes out.</p>	<p>► Check if the appliance is connected to the household outlet. → Disconnect the USB power adaptor from the outlet to use the appliance.</p>
	<p>► This lamp turns off after charging is complete. (See page 16.)</p>
	<p>► Check if the USB power adaptor is fully inserted into the household outlet. Check if the USB plugs are fully inserted into the appliance socket and the USB port.</p>

Problem	Action
<p>My appliance powered off automatically after flashing all the mode indicator lamps (twice per second).</p>	<p>► This indicates that an error shutdown was triggered by the appliance's internal protection feature. Turn the appliance power back on after a while. → If the condition repeats, contact the store where you purchased the unit.</p>
<p>The charging indication/notification lamp () flashes in red (twice per second) when the appliance is being recharged.</p>	<p>► Check if you are recharging the appliance in the recommended recharge temperature range (15 °C to 35 °C). This condition indicates a recharging error. → If the condition repeats, contact the store where you purchased the unit.</p>

Problem	Action
The power turns off during use.	► Charge the appliance. (See page 16.)
	► The power will be automatically turned off after each mode is completed.
	► Check if your hand is in good contact with the Ion inductive panel when applying the head to the skin. (See page 21.)
	► If you switch to a different treatment mode during use, the appliance will turn off automatically approximately 6.5 minutes after the first power-on.

Problem	Action
<ul style="list-style-type: none"> • The appliance cannot be recharged. • The appliance only lasts for a few minutes even after a full recharge. 	► Check if you are using a specified USB power adaptor or USB cable. → Always use the specified power adaptor and cable.
	The battery has reached the end of its life. (approximately 3 years) (See page 31.)
	► Please contact the store where purchased. (Customers cannot replace the battery as doing so may cause malfunction.)
	► Check if the appliance socket is wet. → A wet appliance socket can trigger the USB cable's overheat protection feature to suspend the recharge action. Allow the appliance socket to dry before recharging.

Problem	Action
I cannot feel EMS stimulation.	<p>▶ Refer to “Tips to use EMS effectively”. (See pages 22 to 23.)</p> <p>With the included gel or skin care product you are using applied to the inside of your arm, check to see if you can feel EMS there when the appliance is used in the RF×EMS LIFT mode. (Start checking from level 1.)</p> <p>▶</p>

Battery life


The life of the battery is 3 years. The battery in this appliance is not intended to be replaced by consumers. Have the battery replaced by an authorized service center.

Removing the built-in rechargeable battery

Remove the built-in rechargeable battery before disposing of the appliance.

The battery is to be disposed of safely. Please make sure that the battery is disposed of at an officially designated location if there is one. This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the seal from the bottom of the appliance's main body.
Remove the screw (1 piece) and take out the housing cover ①.
- Remove the screws (2 pieces), take out ② and ③.
Unhook the 3 hooks, remove ③ from ②.
- Insert a ruler between ④, remove the fixture, and take out the battery.

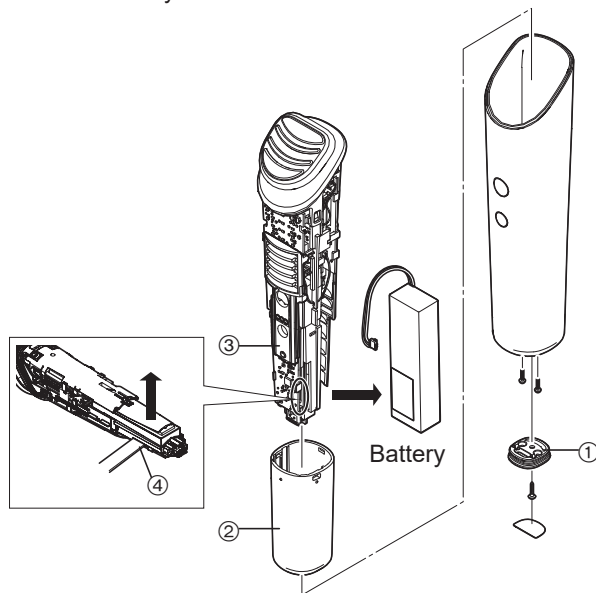For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery. Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Output rating of applicable adaptor (sold separately)	5.0 V $\overline{=}$ 1.0 A In compliance with the Safety Extra Low Voltage (SELV) system
Charging time	Approx. 3 hours

This product is intended for household use only.

Optional parts (Sold separately)

Conductive Gel (EH-4R03)

As this gel does not have to be washed off the skin after use, you can proceed to your regular skin care routine immediately after using the product. Condition the skin with your preferred lotion before applying the gel.

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© Panasonic Corporation 2025
Panasonic Corporation