Panasonic®

Operating Instructions 使用說明書

(Household) Air Massager (家用) 空氣按摩器

Model No. **EW-RAH6**
型號

Leg Refre

English	EN2
繁體中文	TC2

Panasonic Corporation
© Panasonic Corporation 2017
<http://www.panasonic.com>

No.2 EN, TC

Printed in China
中國印刷

W9030RAH422 F1012Y1

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

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What is the Air Massager?

The Air Massager is the household massager that uses airbags for pressure massage. Wrap your legs with the attachments and turn on the manipulator. The airbags swell and massage your legs.


Program	Schematic diagram of massage actions
When your legs feel heavy	Squeeze Powerfully from bottom to top
When you want to relieve your muscle tension by comfortable massage	Knead* Rhythmically
When you want to feel relaxed before going to bed or you feel tired after exercise such as walking	Soft Knead* Softly

• The schematic diagram of the program or massage actions is the representative example.

Please read the contents on Page EN14-EN15 for details.

* Pressure function like knead

Warmth

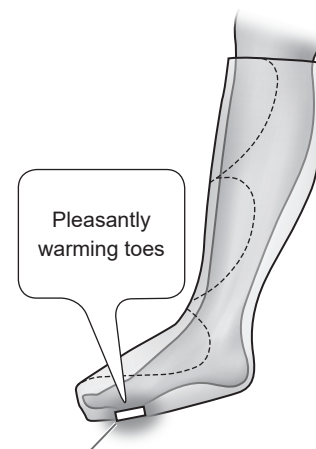Warmer (heater)

• The warmer (heater) warms the textile around toes only.

Contents

Safety Precautions

EN4

Name of each component

EN8

Preparations

EN10

- Wrap
- Use the sole pads / shin pads

Massage

EN14

- To change the program
- Adjust temperature around toes
- To change the intensity

After use

EN16

- Storage
- Keeping

Cleaning and maintenance

EN17

Frequently asked questions

EN18

Before entrusted repair

EN20

Specifications

EN22

Safety Precautions

Always follow these instructions

This section contains instructions that you must follow to prevent personal injury or damage to property.

■ **The instructions are classified according to the level of injury or damage that may occur due to the mistaken use of this unit. Please ensure that you read all instructions before using the Air Massager.**

 **WARNING** Denotes a potential hazard that could result in serious injury or death.

 **CAUTION** Denotes a hazard that could result in minor injury or property damage.

WARNING

Symptoms

● **The following persons should not use this unit.**

- (1) Persons who are prohibited from receiving massages by a physician (e.g.: due to thrombosis, severe aneurism, acute varicose veins, any type of dermatitis and skin infections (including inflammation of the hypodermis), etc.)
- (2) Persons suffering from osteoporosis, spine fracture, sprain or acute pain such as a pulled or torn muscle
- (3) Persons who have a pacemaker, or other personal medical devices prone to electronic interference
- (4) Pregnant women and women who have just given birth
- (5) Persons suffering from severe circulatory disorders in the legs (Otherwise it may cause an accident, ill health or worsening of symptoms.)

● **This unit is not intended for use by persons (including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the unit by a person responsible for their safety. Children should be supervised to ensure that they do not play with the unit.**

(Doing so may cause an accident or injury.)

● **Persons with any of the conditions listed below or who are undergoing medical treatment should consult with their physicians before using the unit.**

- (1) Persons who have a malignant tumor
- (2) Persons suffering from heart disease
- (3) Persons who have a loss of heat sensation
- (4) Persons who have sensory abnormalities due to a severe peripheral circulatory disturbance such as diabetes mellitus
- (5) Persons who have a wound where the massager operates
- (6) Persons whose body temperature is over 38 °C (who have a fever) (E.g.: When a person has acute inflammation symptoms (fatigued, chills, blood pressure fluctuations, etc.) or when debilitated)
- (7) Persons who require bed rest or who are in poor physical shape
- (8) Persons other than those listed above who feel unwell

● **If you begin to feel unwell during use or if you do not feel effects of the air massager, stop use immediately and consult a physician.**

■ **Instructions that you must follow are indicated using the following symbols.**
(The following symbols are examples.)

Denotes a specific operating procedure that must not be performed.

Denotes a specific operating procedure that must be performed.

● **The unit has a heated surface. Persons insensitive to heat must be careful when using the unit.**

(Failure to do so may cause an accident, injury or become unwell.)

When in use

● **Do not use the right attachment for left, and vice versa.**
(Otherwise it may lead to poor body conditions.)

● **Do not use sole or shin pads for parts other than those specified or do not apply them to the bone.**
(Otherwise it may lead to injuries.)

● **Do not stand or walk when in use.**

● **Do not use the unit only for one leg.**

● **Don't use for purposes other than massaging around the calves and feet.**

● **Do not allow children to use the unit, make sure they do not play with the manipulator or attachment.** (Doing so may cause an accident or injury.)

● **Do not use the unit in bedding or kotatsu.**
(Otherwise it may lead to burn or electric shock.)

● **The service time shall be less than 10 minutes per day.**

(Otherwise it may lead to a low-temperature burn*, the opposite effect or injuries.)

* Burns may occur at relatively low temperatures (40 °C to 60 °C) when in contact with the same area of skin for long periods of time, even if you do not feel heat or pain.

The elderly, persons having a delicate skin or insensitive cutaneous sensation may have redness on their skin just after use. If the redness does not disappear after one hour, consult the doctor.

● **Before use, make sure that the cloth is not damaged. Stop using, unplug the AC adaptor and entrust the repairs to qualified persons even for minimum damages.**

(Otherwise it may lead to injuries or electric shock.)

● **For initial use, select "Intensity: Low" mode until you are accustomed to it.** (Otherwise it may lead to injuries or poor body conditions.)

Safety Precautions (Continued)

⚠ WARNING

Power plug, power cord etc.

- **Do not do anything that may damage the power cord or power plug.**
Do not scratch, modify, use near a heater, excessively bend, twist, pull, place heavy objects on top of, or bundle the cord or plug.
(Doing so may cause electric shock, burning or fire due to a short circuit.)
→ Contact the nearest authorized service center for repairs to the power cord or power plug.
- **When you unplug the AC adaptor, do not pull the power cord.**
(Doing so may cause electric shock, burning or fire due to a short circuit.)
- **Do not plug or unplug the AC adaptor or air tube plug with wet hands.**
(Doing so may cause electric shock.)

- **Always use sockets and wiring devices at the correct rated value.**
- **Use 220 V – 240 V ~ (Cannot be used overseas, nor with a transformer.)**
- **Always insert the power plug fully into the socket.**
(Failure to do so may result in electric shock or fire due to overheating.)

If a malfunction or breakdown occurs


- **If a malfunction or breakdown occurs, stop using the unit immediately and unplug the power plug from the wall socket.**
(Failure to do so may result in smoke / fire or electric shock.)

<Examples of malfunctions / breakdowns>

- AC adaptor or air tube plug is damaged
- The plug of the AC adaptor is loose
- The air tube plug is loose when it is inserted to the manipulator

→ Please do not repair it yourself. Entrust the authorized service center with inspection and repairs immediately.

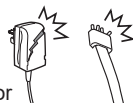Take note of the following points as well

- **Do not stick a pin or needle into the attachment or scar with a knife.**
(Otherwise it may cause fire, an electric shock or burn due to damaged warmer (heater) on the attachment.)
- **Do not modify or disassemble.** (Doing so may cause fire, or may cause the unit to function abnormally resulting in injury.)

⚠ CAUTION

Symptoms

- **If symptoms such as rashes, reddening, itching, etc. occur as a result of using this unit, stop using the unit and consult a physician.**
(Failure to do so may cause an accident or adversely affect health.)

Power supply

- **During cleaning and maintenance, unplug the AC adaptor from the socket.** (Failure to do so may cause an electric shock or burn.)
- **Unplug the AC adaptor from the socket when not in use.** (Since dust and humidity deteriorates insulation, which may result in fire due to electrical fault.)
- **The unit is only to be used with the AC adaptor provided with unit.**


Before and during use

- **Don't use in places with high humidity, such as bathrooms.**
(Doing so may cause electric shock.)
- **Do not immerse manipulator or attachment in water or wash with water.** (Doing so may cause electric shock, burning or fire due to a short circuit.)
- **Do not fall asleep while using the unit.**
- **Do not use the unit after drinking alcohol.** (Doing so may cause an accident or injury.)
- **Don't allow pins, garbage or water etc. to enter the AC adaptor's power plug or the gap of the appliance plug.**
- **When using the massager, don't unplug the AC adaptor from the socket or pull out the appliance plug from the manipulator side.**
(Doing so may cause electric shock, burning or fire due to a short circuit.)
- **When pulling out the air tube plug, do not pull the air tube for connection.**
- **For storage, do not wrap the power cord around the manipulator or the power plug.**
- **Do not fold the air tube for connection or wrap it on the attachment.** (Otherwise, the power cord or air tube for connection may be disconnected due to an overload of the power cord, causing fire or electric shock due to short circuit.)
- **Do not use at the same time as another medical device.**
(Doing so may cause an accident or deteriorated health.)
- **Don't use accessories other than the provided sole and shin pads.** (Otherwise it may lead to injuries.)
- **Do not use AC adaptors other than those specified (Model No.: EW-3A08) or use the provided adaptor for other equipment.**
(Doing so may cause electric shock, burning or fire due to a short circuit.)


- **Unplug the AC adaptor immediately in case of a power failure.**
(Otherwise it may lead to accidents.)

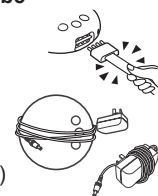
Name of each component

Manipulator

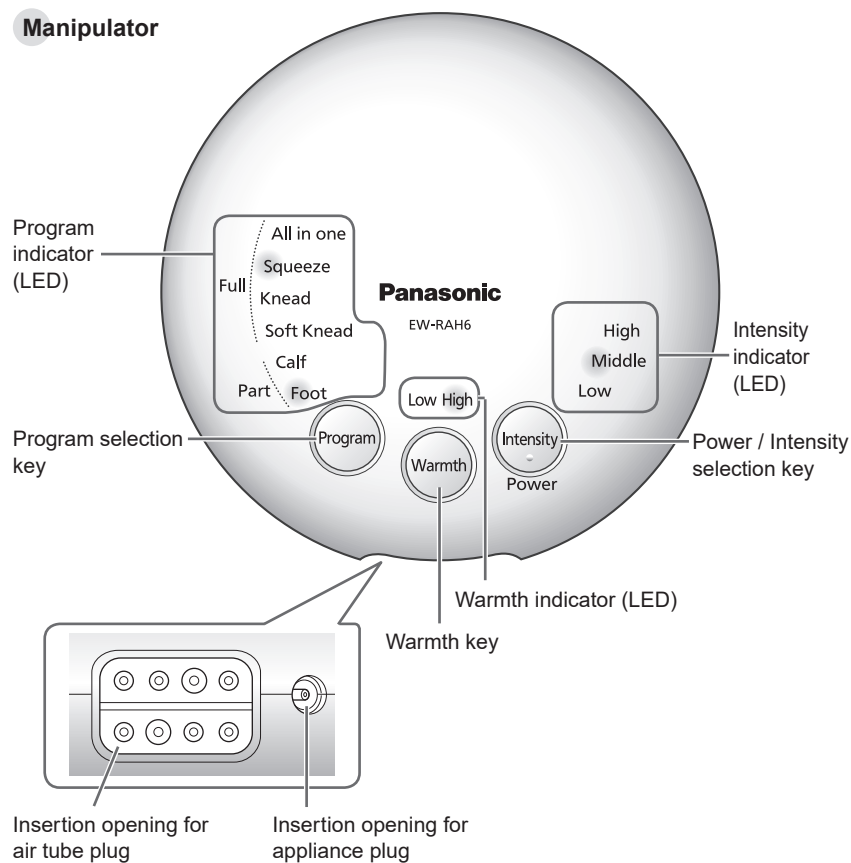AC adaptor

(Model No.: EW-3A08)

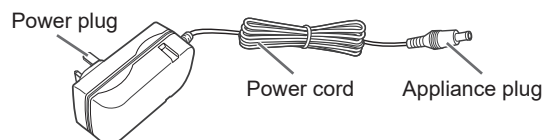
The following symbol indicates that a specific detachable AC adaptor is required for connecting the electrical appliance to the supply mains. The type reference of AC adaptor is marked near the symbol.

Attachment

The unit is composed of right and left attachments. (Bilateral symmetry)


The photo shows the left attachment.

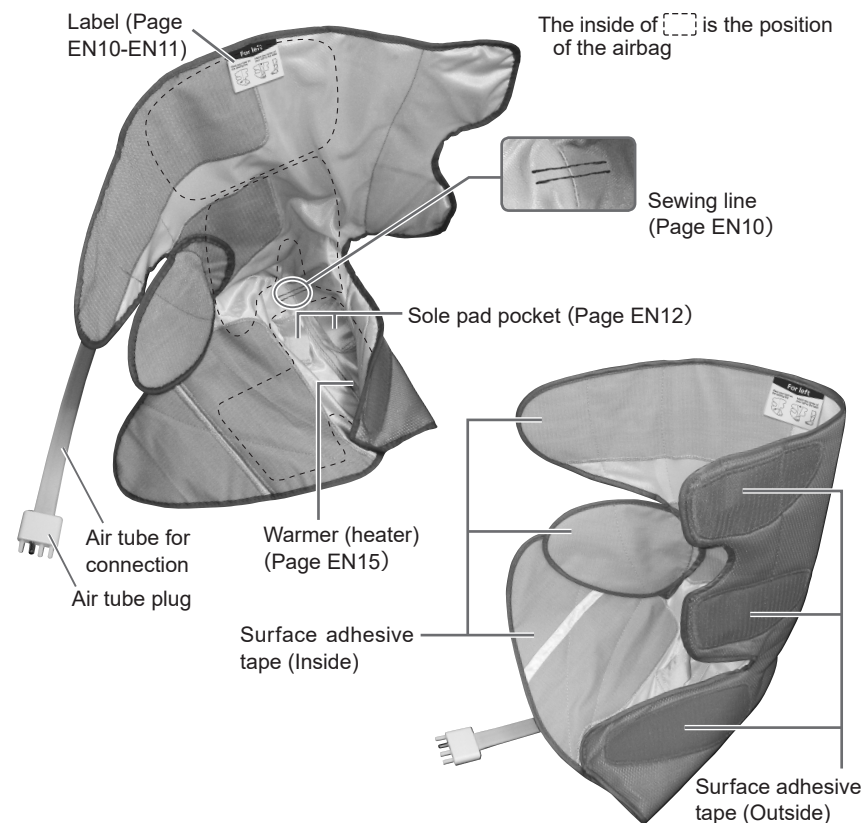Accessories


Sole pads

(2 pcs)

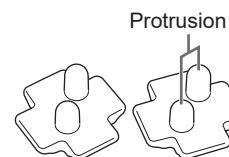Shin pads

(2 pcs)

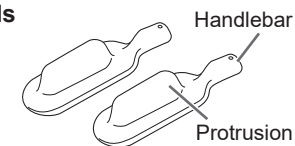Name of each component

Preparations

Wrap

First of all, confirm the following contents

- Whether the cloth is damaged?
- Whether it operates normally after switching on? (Page EN14) (The confirmation is also necessary if it is not used for a certain period of time)
- In case of abnormalities, please refer to Page EN20-EN21.

• **The unit may leave marks after use. If you are concerned about marks...**

Put the unit on over socks and trousers.

The unit leaves marks on calves or feet after use and they may not disappear immediately. If you use the unit before going out, note this point.


• **Be careful not to have the surface adhesive tape contact with your clothes or stockings when putting on the unit.**

• **If you want a stronger massage...**

We recommend that you use sole and shin pads (Page EN12-EN13)

1 Connect the AC adaptor

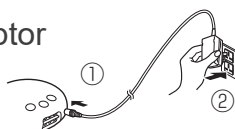
220 V – 240 V ~


2 Check the right and left of the attachments with the label

3 Place your heel on the sewing line

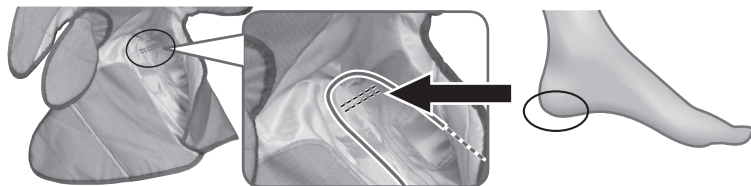
4 Match the center of your calf to the label

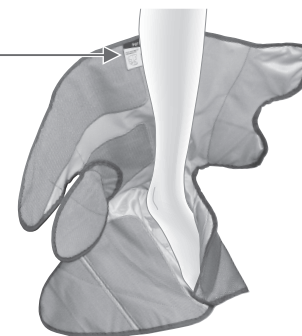5 Wrap the attachments around your legs


- 1) Hold the end of the surface adhesive tape and fix it in order of your **Toe** → **Ankle** → **Calf**
- 2) Press the end of the surface adhesive tape again firmly so that it does not loosen
- 3) Wrap the other attachment around your leg in the same way

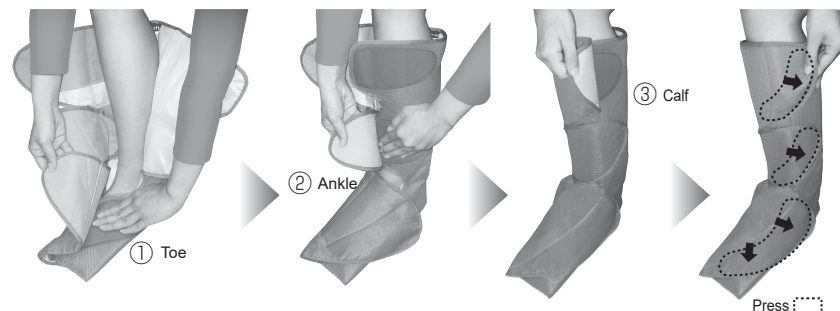- You can adjust the intensity of massage by how you wrap the attachments around your legs.
Wrap them loosely → soft massage
Wrap them tightly → strong massage
- If you put on the unit over thick trousers or leg warmers, the surface adhesive tape may be undone.

6 Insert the air tube plug into the manipulator

After the preparations are completed...

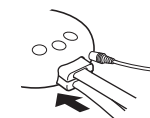Massage (Page EN14)

Preparations



Use the sole pads / shin pads

If you want a stronger massage When you are tired after exercise or walking with high-heeled shoes...

Sole pads


1 Preparation before setting
(Procedure 1 to 2 on Page EN10)

2 Set the sole pads to the pockets

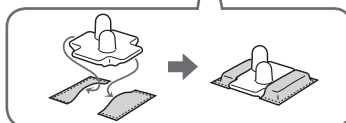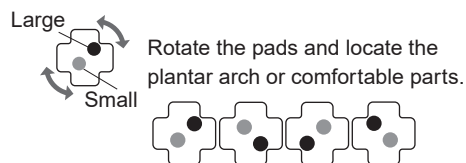Plantar arch

Adjust the large protrusion inside to suit your sole

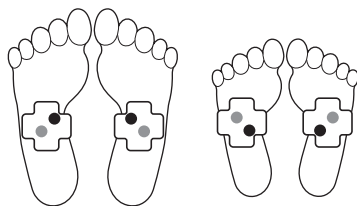Other (comfortable parts)

Locate a comfortable part and adjust the protrusion to it

3 Wrap the attachments around your legs
(Procedure 3 to 5 on Page EN10-EN11)

4 Insert the air tube plug into the manipulator
(Procedure 6 on Page EN11)

- When you use the sole pads, you may be less likely to feel warmth.
- The height of the protrusions are represented with large and small.


Shin pads

1 Preparation before setting
(Procedure 1 to 4 on Page EN10-EN11)

2 Locate the most comfortable position after pressing

- Press around the shin muscles with the shin pad and locate the most comfortable position.

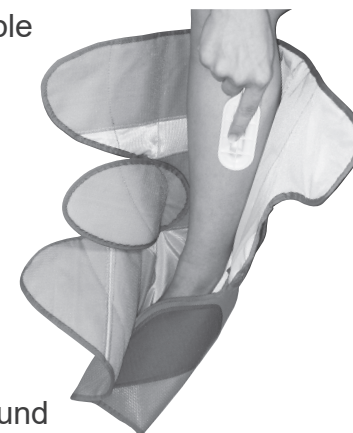3 Wrap the attachments around your legs (Procedure 5 on Page EN11)

4 Insert the shin pad straightly into the part you locate in procedure 2

- Insert the protrusion toward the muscles
- Be careful not to touch the bone
- Set the right and left shin pads in the same way

5 Insert the air tube plug into the manipulator
(Procedure 6 on Page EN11)

Preparations (Use the sole pads / shin pads)

Massage

•To change the program

•Adjust temperature around toes

•To change the intensity

Standard massaging time The service time shall be less than 10 minutes per day

- After using for a long time, it will produce excessive stimulation against the muscles or nerves and sometimes it will lead to the opposite effect.


By pressing , it will switch on automatically

(The unit starts with "All in one", "Warmth: Low", "Intensity: Low")

- After all LED indicators are illuminated in turn, the LED indicator of "All in one", "Warmth: Low", "Intensity: Low" goes on.
- The warmer (heater) begins to warm after 5 minutes.
(The warmer (heater) cannot operate independently without massaging)
- The power will be cut off after about 10 minutes.
(The power will be cut off after about 10 minutes from the start, even when the program is switched over halfway)
- If you want to stop the massage halfway, press  repetitively until it turns to "Off"

WARNING

- For initial use, select "Intensity: Low" mode until you are accustomed to it.
(Otherwise it may lead to injuries or poor body conditions.)

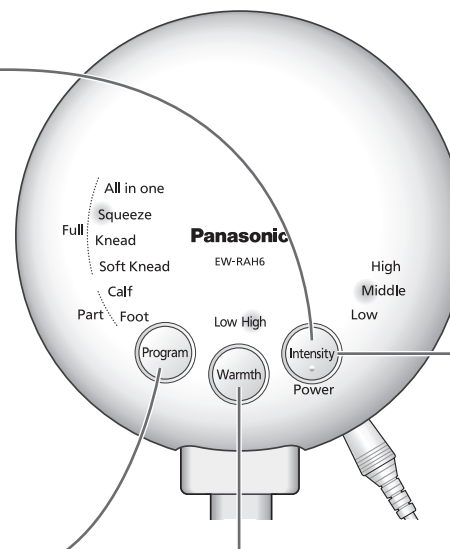To change the intensity


- By pressing the key, the Intensity indicator (LED) will switch in the following order
Low → Middle → High → Off

To change the program

- By pressing the key, the Program indicator (LED) will switch in the following order
All in one → Squeeze → Knead → Soft Knead → Calf → Foot

Adjust temperature around toes

- By pressing the key, the Warmth indicator (LED) will switch in the following order
Low → High → Off

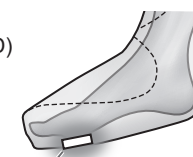Warmer (heater)

- The warmer (heater) warms the textile around toes only.

Program content

Full				Part	
All in one	Squeeze	Knead*	Soft Knead*	Calf	Foot
Specially selected massage with various techniques such as "Squeeze" and "Knead".	Compress the toes to calves like squeezing. Strong massage.	Knead the whole legs fully from the top to the bottom in turn.	Compress the whole of each leg softly and knead them slowly and comfortably.	Complex massage with "Squeeze" and "Knead" focusing on the calves.	Complex massage with "Squeeze" and "Knead" focusing on the toes and soles.

*Pressure function like knead

Massage (To change the program) (Adjust temperature around toes) (To change the intensity)

After use


Be sure that the power supply is in the "Off" state.

Remove the attachments from your legs and unplug the AC adaptor / power plug / air tube plug


Operate according to the steps opposite to those on Page EN10-EN11

Storage

- 1 Fasten the surface adhesive tape

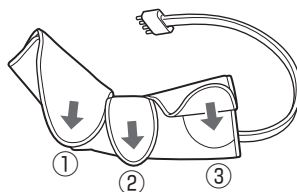
- 2 Bend the air tube for connection lightly and put it in the attachment

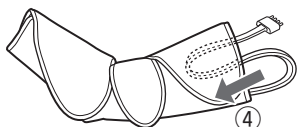CAUTION

- When pulling out the air tube plug, do not pull the air tube for connection.
- For storage, do not wrap the power cord around the manipulator or the power plug.
- Do not fold the air tube for connection or wrap it on the attachment.
(Otherwise, the power cord or air tube for connection may be disconnected due to an overload of the power cord, causing fire or electric shock due to short circuit.)

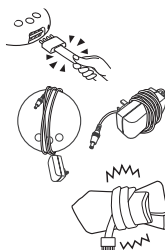Keeping

- Do not place it within the reach of children
- Do not place it at a location with high temperature, high humidity or direct sunlight
- Do not place any objects on the manipulator and attachment surface
(Otherwise the traces may remain on the surface or lead to cross-color or discoloration.)

Cleaning and maintenance


Be sure that the power supply is in the "Off" state.

Do not clean the unit by thinner, benzine or alcohol.
(Otherwise it may lead to malfunctioning, component damage or discoloration.)

- 1 Dip the cloth in water and wring out the water
(For serious dirt, dilute a neutral cleaning agent with warm water or a water solution of about 5%.)

- 2 Wipe the dirt

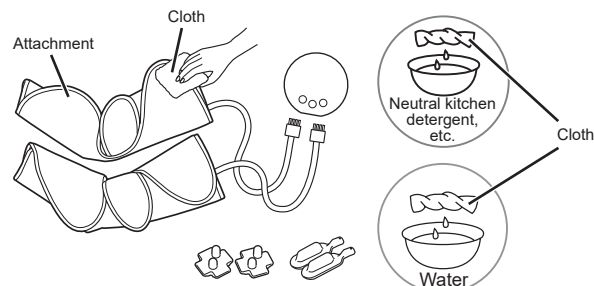- Please avoid using excessive force
- Do not immerse manipulator or attachment in water or wash with water.
(Doing so may cause electric shock, burning or fire due to a short circuit.)

If the surface adhesive tape is stained with dirt or hair, clean it with a toothpick or other tools

- Do not put foreign objects into the insertion opening for air tube plug

Frequently asked questions

Under what circumstances can it not be used?

- Pregnant women and women who have just given birth
- Persons suffering from severe circulatory disorders in the legs

You cannot use the unit in the above cases. (Page EN4)

If you have a fever or symptom in your body, consult doctor before use.

Can I use it many times a day?

- It can be used once in a day. (The power will be cut off after about 10 minutes)
If the time for a single use is less than 10 minutes, the total time for multiple uses shall be less than 10 minutes per day.
After using it for a long time, it will produce excessive stimulation against the muscles and nerves and sometimes it will lead to the opposite effect.

WARNING

- **The service time shall be less than 10 minutes per day**
(Otherwise it may lead to a low-temperature burn, the opposite effect or injuries.)

Can it be used for parts other than calves and feet?

- It can only be used for around the calves and feet.
If it is used for other parts, it may lead to accidents or injuries.

What is the size around the foot and around the calf for the unit?

- Available foot size: approx. 20 cm ~ 28 cm
- Available calf girth: approx. 27 cm ~ 48 cm

The toes are too hot...

- Set the warmth to "Low". If it is still hot, put on thick socks.

The ankles and calves are not warm...

- The warmer (heater) warms the textile around toes only.

The massage intensity is too high...

The massage intensity is too low...

- The massage intensity can be adjusted by the wrapping method of the attachment.
Too high intensity → wrap the attachment loosely.
Too low intensity → wrap the attachment tightly.
- Try to use the sole and shin pads when the massage intensity for soles and shins is weak. (Page EN12-EN13)

I feel uncomfortable when the skin is in contact with the surface adhesive tape...

- Cover the skin with a thin handkerchief or other article.

The air in the attachment cannot be discharged completely...

- If there is a lot of air remaining in the attachment (airbag), it will take some time to discharge. Remove the air tube plug and wait a moment. The air can be discharged on its own.

It sends out the following noise. Is it abnormal?

- Crackling sound
(Tensioning sound of surface adhesive tape when the air enters into the attachment (airbag).)
- Kachi sound
(It is the sound sent out when the manipulator adjusts the air.)
- Humming sound
(Motor sound.)

It is the sound sent by the unit structure. This is normal.

Can it be used abroad?

The unit adopts the power supply of 220 V – 240 V ~.

In addition, when using a transformer or a conversion plug, a malfunction may occur.

Do not use them.

Before entrusted repair

Please confirm the following.

After confirmation, if the abnormality remains, first unplug the AC adaptor from the socket and then entrust to an authorized service center for inspection and repairs.

Symptoms	Items to be confirmed
The power cannot be switched on	<ul style="list-style-type: none"> Does the AC adaptor fall off? (Page EN10)
It stops halfway	<ul style="list-style-type: none"> Whether the Power / Intensity selection key is in the "Off" state? (Page EN14-EN15) Does the air tube plug fall off? (Page EN11)
The air cannot enter	<ul style="list-style-type: none"> Whether the air tube for connection is folded?
Toes do not feel warmth	<ul style="list-style-type: none"> The warmer (heater) warms in 5 minutes after turning on the unit. Is the Warmth indicator (LED) on? Pressing the Warmth key lights the Warmth indicator (LED). (Page EN15) Are your heels adjusted to the specified position? (Page EN10) Are you using the sole pads? (Page EN12) Are you putting the unit on thick socks?
I feel warmth differently on right and left	<ul style="list-style-type: none"> Does one air tube plug fall off? (Page EN11) Are your heels adjusted to the specified position? (Page EN10)

Symptoms	Items to be confirmed
The massage intensity is low	<ul style="list-style-type: none"> Is the attachment wrapped loosely? (Page EN11) Are you putting the unit on thick socks? Does the air tube plug fall off? (Page EN11) Is the air tube plug inserted completely? (Page EN11) Whether the air tube for connection is folded?
The massage intensity is high	<ul style="list-style-type: none"> Is the attachment wrapped tightly? (Page EN11)
Intensity indicator (LED), Program indicator (LED) or Warmth indicator (LED) does not go on	<ul style="list-style-type: none"> Does the AC adaptor fall off? (Page EN10) Press the Power / Intensity selection key, Program selection key and Warmth key and confirm. (Page EN14-EN15)
The air cannot be discharged completely	<ul style="list-style-type: none"> Remove the attachment and unplug the air tube plug. (Page EN16)

In case of abnormalities ...

Manipulator or AC adaptor is overheating	<p>These conditions are failures. To prevent an accident, stop using immediately and contact the dealer.</p>
It is damaged	
The attachment cloth is damaged	

Before entrusted repair

Specifications

Product name		Air Massager
Product model		EW-RAH6
Power supply		220 V – 240 V ~ 50 – 60 Hz
Power consumption		21 W (Rated power of the heat transfer product 5 W)
Air pressure		Low : approx. 35 kPa Middle : approx. 42 kPa High : approx. 50 kPa (It slightly varies depending on the program or part)
Duration		Approx. 10 minutes
Dimension (packing qty)	Manipulator (1 pc)	Approx. 145 mm (D) × 50 mm (H)
	Attachment (2 pcs)	Approx. 570 mm (W) × 300 mm (L) × 390 mm (H)
Weight		Approx. 1 600 g (manipulator + attachment)
Length of power cord		Approx. 1.8 m
Corresponding perimeter of attachment		Foot size: approx. 20 cm ~ 28 cm Calf girth: approx. 27 cm ~ 48 cm
Surface temperature of the inside textile of the attachment (in indoor temperature of 20 °C)		Low: approx. 37 °C High: approx. 40 °C
Manufacturer		Panasonic Manufacturing (Shanghai) Co., Ltd. No. 258 East Jiangnan Road, Songjiang Industrial Park, Shanghai
Country of origin		China

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