Panasonic®

Operating Instructions 使用說明書


(Household) High-Frequency Conductive Belt
(家用) 高周波傳導帶

Model No.
型號

EW-RAE50

English EN2

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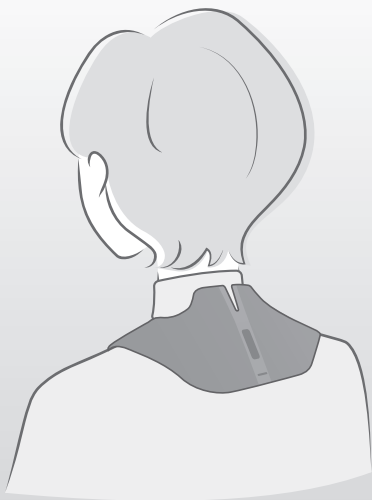Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

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
Household High-Frequency Conductive Belt

High-frequency pulses penetrate deep into body tissues, promoting comfortable shoulder relaxation and alleviating fatigue.

Purpose or Effect of Use

Used to relieve shoulder fatigue. For general home use.

Generates high-frequency pulses from 12 antenna points

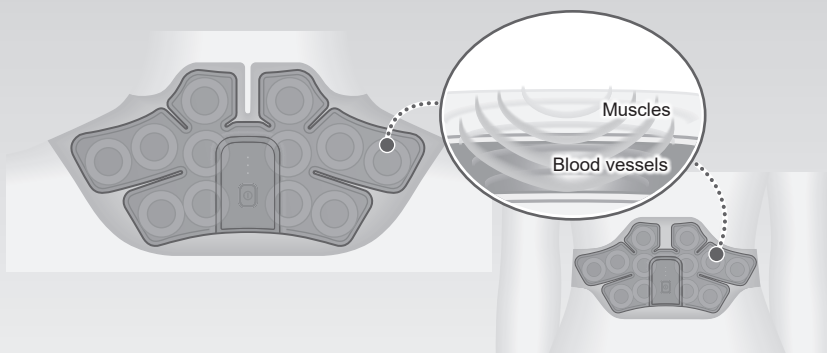Can be worn for extended periods

- Because it uses non-stimulating high-frequency pulses, long-term use is possible.
- Lightweight design; operates continuously for approximately 12 hours.

Difference from Low-Frequency Products

(Comparison with our company's products)

	High Frequency	Low Frequency
Structure	Frequency above 5 MHz Delivers high-frequency electromagnetic energy onto uncomfortable shoulder areas to relieve fatigue	Frequency around 3-1200 Hz Allows the target area to relax via gentle low-frequency currents on the skin surface
Stimulation to the body	No sensation of stimulation	Sensation of stimulation
Features	Non-irritating and suitable for extended use. Continuously relieves fatigue while worn	Short duration (approx. 15 minutes), stimulation is physically felt

Use while working

Use while Doing Housework

The main unit can be placed inside a belt (sold separately) for use on the lower back.

Belt (Model: **EW-9RE5W**)

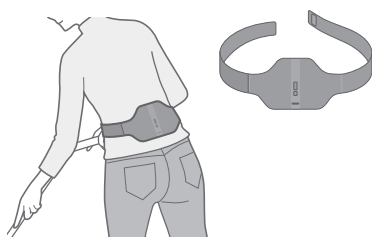
Table of Contents

Important Safety Instructions **EN4**

Part Names and Accessories **EN10**

Charging **EN12**

Confirming High-Frequency Pulse Activation **EN14**

Usage after Wearing **EN15**

After Use **EN16**

Maintenance **EN17**

Storage **EN18**

•Disposal

Belt (Sold Separately) **EN19**

FAQ **EN20**

Before Requesting Repair **EN22**

Consumables **EN23**

Specifications **EN24**

Warranty and After-Sales Service **Back Cover**

Important Safety Instructions

Always follow these instructions

This section contains instructions that you must follow to prevent personal injury or damage to property.

- The instructions are classified according to the level of injury or damage that may occur due to the mistaken use of this unit.

	DANGER	Denotes a potential hazard that could result in serious injury or death.
	WARNING	Denotes a potential hazard that could result in serious injury or death.
	CAUTION	Denotes a potential hazard that could result in minor injury or property damage.

- Instructions that you must follow are indicated using the following symbols.

	Denotes a specific operating procedure that must not be performed.
	Denotes a specific operating procedure that must be performed.

DANGER

Symptoms

-  ● Do not use simultaneously with the following medical electrical equipment.
 - (1) Persons wearing implantable medical electronic devices such as pacemakers that are susceptible to electromagnetic interference.
 - (2) Wearable medical electrical equipment such as electrocardiographs.
 - This appliance shall not be used by, or on, persons wearing implants not compatible with its use.
- (Failure to do so may lead to accidents or significant physical injury.)

WARNING

Symptoms

-  ● Do not wrap around areas other than the shoulders
 - This appliance is not suitable to treat certain parts of the human body.
 - Persons (including children) with physical, sensory or mental disabilities, or lacking operational experience or knowledge, shall not use this unit without supervision and guidance regarding the use of this unit from a guardian. Prevent children from treating this product as a toy, as this may lead to danger.
- (Failure to do so may lead to accidents, injuries, or deterioration of health.)

WARNING

Symptoms

- This appliance shall not be used by or on persons up to and including 14 years of age.
Additionally, do not allow children to play on or sit on the unit.
(Failure to do so may lead to accidents or injuries.)
- Use of this appliance may cause discomfort or adverse effect to, persons suffering from skin disease or inflammation of the skin, persons with sensitive skin and persons having an allergic constitution.
- This appliance is not to be used on the parts where the integrity of the skin is damaged.

- Persons with any of the conditions listed below or who are undergoing medical treatment should consult with their physicians before using this unit.
 - (1) Persons who have an acute disease
 - (2) Persons who have a malignant tumor
 - (3) Persons suffering from heart disease
 - (4) Persons who are during menstruation, pregnancy, or shortly after delivery
 - (5) Persons who have sensory abnormalities due to a severe peripheral circulatory disturbance such as diabetes mellitus
 - (6) Persons who have an infectious disease
 - (7) Persons suffering from anemia
 - (8) Persons suffering from skin sensory impairment or skin abnormalities
 - (9) Persons who have an abnormal blood pressure (≥ 180 mmHg)
 - (10) Persons who have acute conditions such as spinal fractures, sprains, or muscle strains
 - (11) Persons whose body temperature is over 38 °C (who have a fever) (E.g.: When a person has severe acute inflammation symptoms [fatigued, chills, blood pressure fluctuations, etc.] or when debilitated)
 - (12) Persons who require bed rest or who are in poor physical shape
 - (13) Persons other than those listed above who feel unwell
- Discontinue use immediately and consult a physician if you experience physical discomfort or fail to perceive the intended therapeutic effects.
(Failure to do so may lead to accidents or injuries.)
- Children must be strictly supervised to prevent them from playing with this unit.
(Failure to do so may lead to accidents or injuries.)

Important Safety Instructions (Continued)

WARNING

Before and during use

- Do not use simultaneously with other therapeutic devices.
Also, do not use with topical medications (including sprays, etc.)
(Failure to do so may lead to accidents, injuries, or deterioration of health.)
- Do not use in high-humidity areas such as bathrooms.
(Failure to do so may lead to electric shock accidents or fire caused by short circuits.)
- Do not use while sleeping.
(Failure to do so may lead to accidents or deterioration of health.)
- Do not bypass or remove the inherent safety device of the appliance.

- Before use, be sure to check the unit for damage.
If damaged or not working properly, stop using immediately and send the unit for repair.
(Failure to do so may lead to accidents or deterioration of health.)

If a malfunction or breakdown occurs

- Stop using immediately and unplug the power plug if an abnormality or malfunction occurs.
(Failure to do so may lead to smoke, overheating, fire, or electric shock.)

<Examples of malfunctions/breakdowns>

Main unit, USB power adapter, dedicated USB cable

- Unit is damaged, deformed, or abnormally hot.
- Emits a burning smell.
- Power plug is not fully inserted.
- USB cable is not fully inserted.

→ Do not attempt to repair it yourself; immediately contact an authorized service center for inspection and repair.

WARNING

Please observe the following points

- Do not wash with water or drip water onto the main unit (except for the strap).
(Failure to do so may lead to overheating, rupture, or fire caused by short circuits.)

- Do not modify or disassemble the unit without authorization.
(Failure to do so may lead to accidents, injuries, or health deterioration caused by overheating, rupture, fire, electric shock, leakage, or abnormal operation.)
- Do not charge, use, leave, or store in high-temperature environments (near fire, under strong sunlight, inside a car, etc.).
- Do not place the main unit in fire or heat it.
- Do not pierce the main unit with a nail or subject it to impact.
(Failure to do so may lead to injury or environmental contamination due to overheating, rupture, fire, or leakage of the built-in battery.)
- Do not store within reach of infants or young children.
(Failure to do so may lead to accidents or injuries.)

- This appliance contains batteries that are only replaceable by skilled persons.
- Only use consumable materials recommended by the manufacturer.
- This appliance is only to be used with the unit provided, the USB Power Adapter Model No. **SW-050100B**
(Failure to do so may lead to accidents or injuries.)


Important Safety Instructions (Continued)

WARNING

Regarding USB Power Adapter and Dedicated USB Cable

- Do not use USB power adapters or dedicated USB cables other than those specified, and do not use them for other products.
(Failure to do so may lead to accidents, injuries, overheating, rupture, electric shock, or fire caused by short circuits.)
- Do not plug or unplug the power plug or USB plug with wet hands.
(Failure to do so may lead to electric shock.)

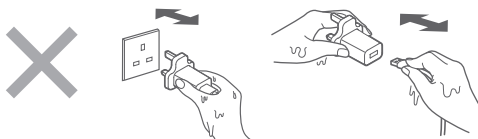- Do not allow pins, debris, water, etc., to enter the gaps of the main unit, signal device, USB power adapter, or USB cable.
(Failure to do so may lead to electric shock or fire caused by short circuits.)
- Do not perform any action that may damage the USB power adapter or dedicated USB cable.
Do not scratch, modify, place near heaters, excessively bend, twist, pull, pinch, place heavy objects on, or bundle them.
(Failure to do so may lead to electric shock or fire caused by short circuits.)
- Consult the store where you purchased the unit or an authorized service center regarding how to repair the power adapter.
- Do not pull on the cord when unplugging the dedicated USB cable.
(Failure to do so may lead to electric shock or fire caused by short circuits.)

- Be sure to insert the power plug and USB plug all the way in to ensure they are securely seated.
(Failure to do so may lead to electric shock or fire caused by overheating.)

Power Supply

- Use sockets and wiring devices with correct rated values.
- Use 220 V AC. (Do not use this product abroad, and transformer usage is prohibited.)
(Failure to do so may lead to electric shock or fire due to overheating.)

CAUTION

Before and during use

-  ● Do not use tape to stick the main unit directly onto the skin.
(Failure to do so may lead to skin inflammation.)
- Do not press firmly against the body or use while rubbing against something.
(Failure to do so may lead to internal bleeding, abrasions, etc.)
-  ● If rash, redness or itching occurs during use of this unit, discontinue use of this unit and consult a physician.
(Failure to do so may lead to accidents or physical discomfort.)
- Limit continuous use to a maximum of 12 hours.
(Failure to do so may lead to skin inflammation or physical discomfort.)

Important Safety Instructions

USB Power Adapter, Dedicated USB Cable, etc.

-  ● Do not wrap the cable around the main unit, signal device, or USB power adapter, or twist it.
(Failure to do so may strain the cable, leading to wire breakage, resulting in electric shock or short circuit.)


-  ● When cleaning and maintaining, be sure to unplug the USB plug and the power plug.
(Failure to do so may lead to electric shock, injury, or overheating, rupture, fire, leakage of the built-in battery causing injury or environmental contamination.)
- When not in use, be sure to unplug the power plug from the outlet.
(Failure to do so may lead to insulation deterioration due to dust or moisture, leading to electrical leakage or fire.)
- In the event of a power outage, unplug the power plug immediately.
(Failure to do so may lead to accidents or injuries when power is restored.)


Part Names and Accessories

When taking the unit out of the color box, the main unit is already placed inside the shoulder strap.


Main Unit

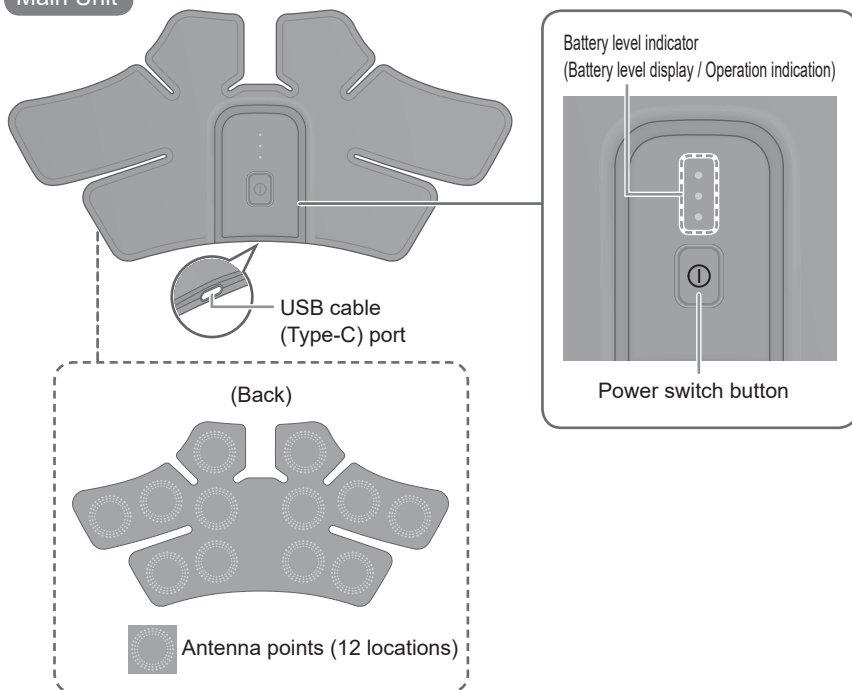
USB Power Adapter (Accessory) (Model: SW-050100B)

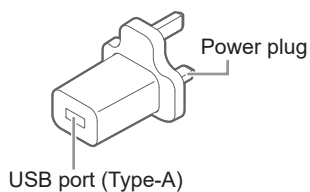
Dedicated USB Cable (Accessory) (Model: EWRA550W7697)

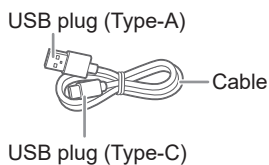Shoulder Strap (Accessory)

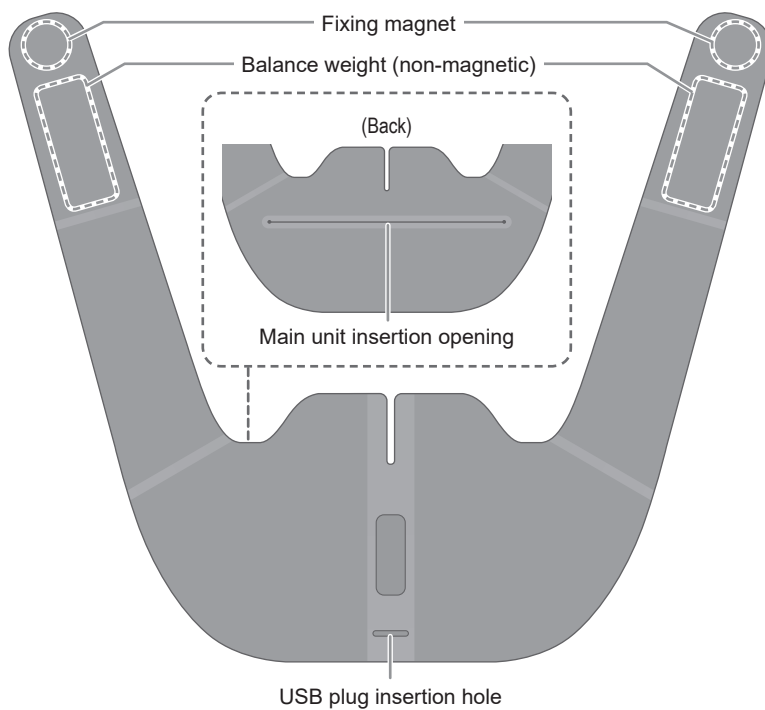Signal Device (Accessory)

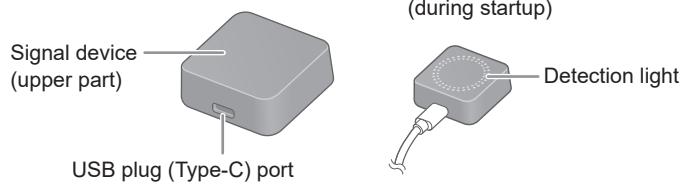
Charging

The battery is not fully charged at the time of purchase. Take out the main unit from the shoulder strap and be sure to charge it before use.


1 Removing the Main Unit from the Shoulder Strap

- ① Remove the main unit from the insertion opening on the back of the shoulder strap.

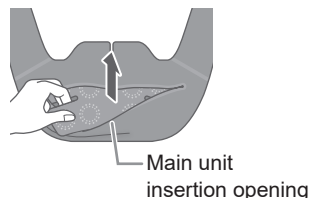2 Charging

- ① Insert the dedicated USB cable into the main unit.
- ② Insert the power plug into the outlet

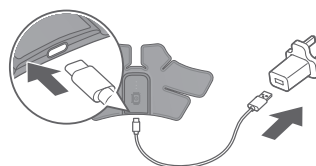Indicator Display During Charging (On Flashing Off)

Charging (Battery level below 25%)	Charging	Charging	Charging Complete (After 3 seconds)


- Charging completes in approx. 2 hours (charging time depends on remaining battery level).
- Continuous usage time when fully charged: Approx. 12 hours
- Charging time varies depending on remaining battery level. Charging for more than 2 hours does not affect the main unit.

3 After the battery level indicator turns off, unplug the power plug.


- Unplug in the order of ② then ① as shown in Step 2.

4 Inserting the Main Unit into the Shoulder Strap


- ① Align the back of the shoulder strap with the back of the main unit, matching top and bottom

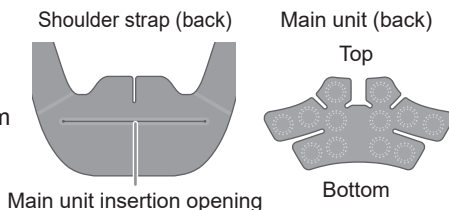- ② Open the main unit insertion opening and insert the main unit.

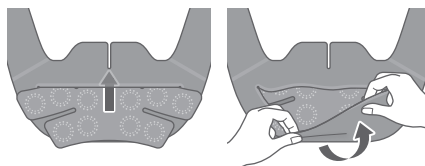- ③ Ensure the front side is facing up and the power switch button on the main unit is visible

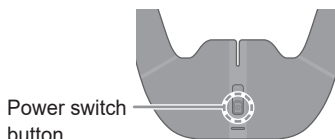Charging

- Maximum output power of the product: 29.52 mW
- Maximum output energy of the product: 106.272 J
- If the usage time is very short even after charging correctly, the built-in battery may have reached the end of its life. Contact an authorized service center for inspection and repair. (The battery life is approximately 3 years, though it may vary depending on actual usage and storage conditions.)
- Recommended charging temperature is 5 °C to 35 °C. Charging outside the recommended temperature range may cause deterioration of the built-in battery or inability to charge.
- If the battery level indicator flashes rapidly (approx. twice per second), it may indicate a malfunction due to charging outside the recommended temperature range. Charge within the recommended charging temperature range.
- Charging time may vary when using for the first time or after not using for more than six months. The indicator may not light up or flash for a few minutes after charging starts, but it will eventually light up or flash if charging continues.

EN13

Confirming High-Frequency Pulse Activation


You can confirm using the signal device. (Via buzzer sound and signal device detection light display)

1 Starting the signal device

- ① Insert the dedicated USB cable into the signal device and USB power adapter.

- ② Insert the power plug into the outlet.

- The signal device detection light turns on for 1 second, then turns off, indicating the signal device has started.

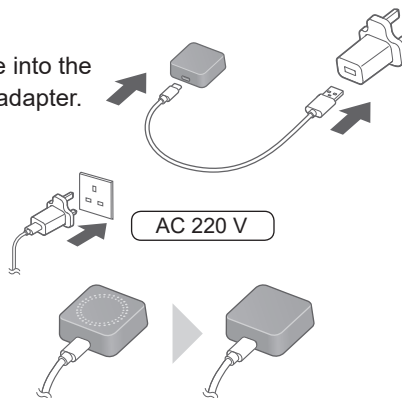
2 With the main unit power turned on Bring the antenna part close to the signal device

- When high-frequency pulses are detected, the signal device emits a buzzer sound. (Buzzer sound and signal device detection light flashing)

<Regarding Reaction>

Distance to Main Unit	Buzzer Sound	Signal Device Detection Light
Far	Beep, beep, beep...	Flashing slowly
Near	Beep beep beep beep...	Rapid Flashing

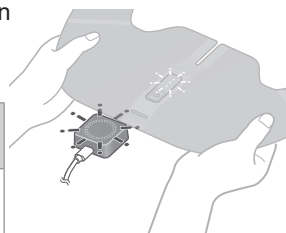- Depending on the positional relationship between the main unit and the signal device, It is easier to trigger a response if the signal detector detects high-frequency pulses from multiple antennas.
- It may be difficult to get a response on a metal table.

3 Unplugging the Power Plug

- Unplug in the order of ② then ① as described in Step 1.

Usage after Wearing

Wear over thin clothing when using.

Please confirm


Does the battery level indicator light up or flash normally after turning on the power?
(Confirm even when not using temporarily.)

→ Please check EN22-EN23 if any abnormality occurs.

1 Short Press the Power Switch Button


- The battery level indicator lights up for 3 seconds, the indicator flashes.

(Display example: When fully charged)

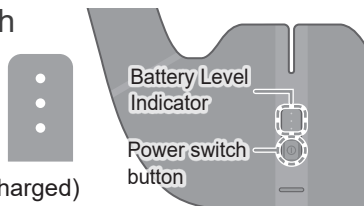
Indicator Display During Use

(Flashing  Off )

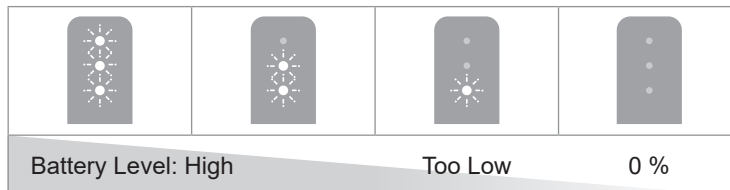- All indicators turn off when the battery has no remaining charge.
(Power off)

2 Position the battery level indicator on the outside of the product, Wear along the neck and shoulders (main unit in the middle of the back)

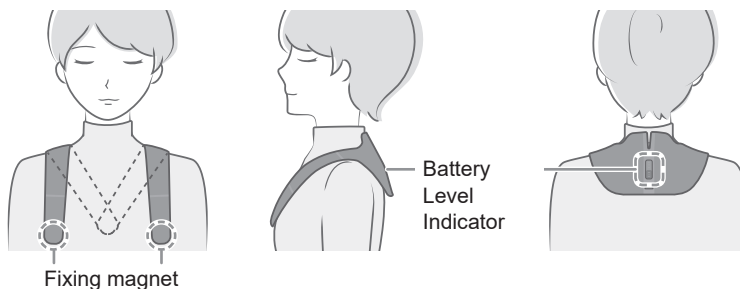
- Use the magnet for fixation if you wish to secure the unit's position or are worried about it falling off while going out.


After Use

- 1 Remove the shoulder strap from the shoulder,
Short press the power switch button to turn off the power

- The battery level indicator turns off.

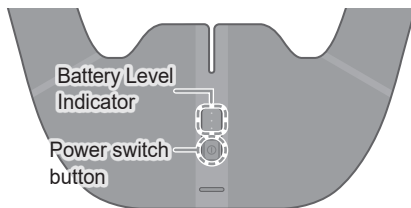- Charging can also be performed while the main unit is inside the shoulder strap.

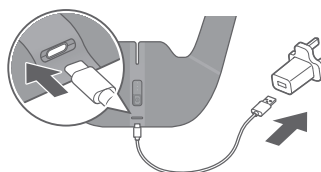Maintenance


Be sure to unplug the USB plug and remove the main unit from the shoulder strap before performing the following operations.

Do not use chemicals such as thinner, benzene, or alcohol
(as this may cause component malfunction, cracking, discoloration, etc.)

1 Wipe with a dry cloth

- Do not wash with water. (Except for the shoulder strap)
- Do not subject the main unit or signal device to excessive force or drop them.
(Failure to do so may lead to malfunction.)

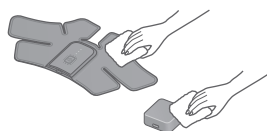After Use

The shoulder strap is fully washable.

Washing Method

- ① Put into a laundry net.
- ② Soak in water below 30°C and wash by gently pressing.
 - Do not use detergents containing bleach.
(Failure to do so may lead to discoloration.)
 - Do not rub or scrub. (Failure to do so may lead to fabric damage.)
 - Wash separately from other items. (In order to prevent color transfer.)
- ③ Rinse thoroughly.
- ④ Gently press to drain water.
 - Do not wring or twist to dehydrate.
(Failure to do so may lead to fabric damage or deformation.)
- ⑤ Gently smooth out wrinkles and dry in the shade.
 - Do not expose to direct sunlight or use a tumble dryer.
(Failure to do so may lead to fabric shrinkage or discoloration.)
 - Dry immediately after washing. (To prevent color transfer.)
 - Do not iron. (Failure to do so may lead to fabric damage or deformation.)

Composition and Washing Instructions

Nylon	80 %
Spandex	20 %

Maintenance

Storage

Long-term storage may cause deterioration of the built-in battery.
Be sure to fully charge before storage.

- Do not place in the following locations:
 - Locations with high temperature or humidity
(Failure to do so may lead to malfunction or fabric discoloration.)
 - Locations exposed to direct sunlight
(Failure to do so may lead to fabric discoloration.)
 - Locations accessible to infants and young children
(Failure to do so may lead to accidents or injuries.)
- Do not place other items on the main unit or accessory strap.
(Failure to do so may lead to indentations, color transfer, or discoloration on the surface.)
- Do not move while holding the cable when the dedicated USB cable is connected.
(Failure to do so may lead to injury or malfunction.)
- Do not forcibly fold or roll up the main unit. (Failure to do so may lead to malfunction.)
- Keep away from pets such as dogs and cats.
(Failure to do so may lead to damage or malfunction.)

Disposal of used batteries

The built-in battery is a valuable recyclable resource.
Before disposing of this unit, the battery must be removed from the main unit.
The main unit must be powered off when removing the battery. Batteries should be disposed of safely.

Removal steps: 1. Use a professional tool to pry open the silicone operation part. 2. Remove the fixing screws. 3. Remove the battery.
(The used battery will be handed over to the recycling agency set up by the government for unified disposal; Or the battery is handed over to professional institutions and their technicians for recycling.)

The built-in battery has passed CCC mandatory certification.

Disposal instructions for the following consumables:

- USB Power Adapter
 - Shoulder Strap
 - Dedicated USB Cable
 - Signal Device
- (Please dispose of in the recycling bin)

Belt (Sold Separately)


The belt is sold as an accessory.
Purchase separately. (EN23)

Applicable size:

Waist circumference approx. 58-100 cm
(Using extension belt)

Wear over thin clothing when using.



- Refer to EN12-EN16 for details on charging and operation of this unit.

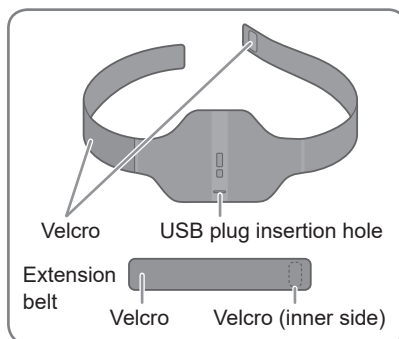
1 Inserting the Main Unit into the Belt

- ① Align the belt with the back of the main unit, matching top and bottom.


Waist accessory belt (back) Main unit (back)
Main unit insertion opening Top Bottom

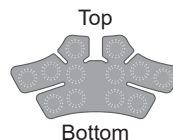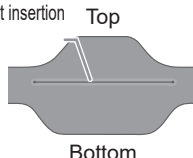- ② Open the main unit insertion opening and insert the main unit.

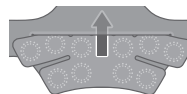
- ③ Ensure the front side is facing up and the power switch button on the main unit is visible.

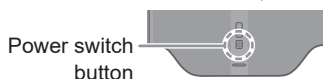2 With the main unit power turned on Wrap around the waist and secure with Velcro

- If the belt length is insufficient, use it in combination with the extension belt.

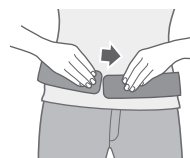3 After use, remove from the waist.

- Turn off the power of main unit after removal.

- The belt is washable. Take out the main unit, fasten the Velcro, and wash according to the instructions on EN17.

- If debris or hair gets tangled in the Velcro, use a toothpick or similar tool to remove it.

Storage

Belt (Sold Separately)

EN19

FAQ

Do you feel nothing?

- The conductive belt provides no stimulation. (EN2)

How long does it take to see effects?

- Effects vary depending on the individual and the severity of pain; continuous long-term use is recommended.

It is difficult to see the flashing of the battery level indicator, and I am unsure if high-frequency pulses are being output.

- Even in bright rooms, it may be difficult to see.
Use the signal device to check if high-frequency pulses are being output. (EN14)

The area where it is worn feels warm; does it have a heating function?

- It may feel warm after wearing, but this unit does not have a heating function.

Can I use it while sleeping?

- Do not use while sleeping.

Can I use a USB power adapter or USB cable other than the supplied ones for charging?

- Be sure to use the USB power adapter and dedicated USB cable supplied with this unit to avoid accidents, injuries, or malfunctions.

What is the waist size for the belt (sold separately)?

- Applicable size: Waist circumference approx. 58-100 cm (using the extension belt) (EN19)

Can the shoulder strap or belt (sold separately) be washed in a washing machine?

- Please wash by hand. (EN17, 19)

Can it be used overseas?

- This product is for use in Hong Kong only and cannot be used overseas. Also, using a transformer or plug adapter may cause malfunction.

Does it affect anti-theft systems?

- The high-frequency pulses from the main unit may trigger a reaction when approaching anti-theft gates in stores or airports. In such cases, turn off the power of the main unit.

Can it be left unused for a long time?

- Since the built-in battery may deteriorate, please fully charge it every six months if you do not use it for more than half a year.

Does frequent charging affect the lifespan of the built-in battery?

- A battery that does not easily deteriorate even with frequent charging is used, so it will not affect the battery life.

FAQ

Before Requesting Repair

Please check “FAQ” (EN20-EN21) and the following contents.

If abnormalities persist after checking, first unplug the USB power adapter from the power outlet and contact the retailer where you purchased the product. (Refer to EN25 for details)

Symptoms	Please confirm
Fail to charge. (Battery level indicator does not light up or flash.)	● Is the USB cable unplugged? (EN12) ● Can I charge using a USB power adapter or dedicated USB cable other than those supplied with this unit? (EN10, 20) ● Charging time may vary when using for the first time or after not using for more than six months. The indicator may not light up or flash for a few minutes after charging starts, but it will eventually light up or flash if charging continues.
All battery level indicators flash rapidly (approx. twice/second)	● Is charging being performed within the recommended charging temperature range (5-35°C)? (EN13)
After charging is complete, the dedicated USB cable is inserted into the main unit to charge again.	● Even if charging is complete, the battery level indicator will light up for the first 3 seconds so you can clearly confirm the connection of the dedicated USB cable to the main unit and the remaining battery level. This is not a malfunction.
Cannot use for a long time even after charging correctly	● Lifespan of the built-in battery. (EN13)
The battery level indicator does not light up or flash when power is turned on / No beep sound when the antenna is close to the signal device	● Is the battery depleted? (EN15) ● Is the main unit's power turned on? (EN15) ● Is the dedicated USB cable unplugged from the signal device? (EN14)
Beeping occurs even when the main unit is not close to the signal device	● It may be affected by surrounding electronic devices; this is not a malfunction. Move the signal device away from electronic equipment.

In such cases...

Symptoms	Please confirm
Main unit, signal device abnormally overheating ▶	Malfunction. To prevent accidents, stop using and contact a Panasonic authorized service center for repair.
Damaged ▶	
Abnormal noise or odor ▶	

Consumables

USB Power Adapter

Model: **SW-050100B**

Shoulder Strap

Belt

Dedicated USB Cable

Signal Device

Before Requesting Repair

Consumables

Consult the store where you purchased the product regarding the signal device.

Specifications

Product Name		High-Frequency Conductive Belt
Product Model		EW-RAE50
Rated Parameters	During Main Unit Operation	Rated Voltage DC 3.8 V
	During Main Unit Charging	Power Supply AC 220 V 50 Hz
		Rated Power Consumption 0.4 W
	During High-Frequency Check Program	Power Supply AC 220 V 50 Hz Rated Power Consumption 0.8 W
Category		General Household Appliances
Common Name		Household High-Frequency Conductive Belt
Executive Standard		IEC60335-1, IEC60335-2-115
Dimensions (D × W × H)	Main Unit	Approx. 236 x 117 x 8 mm
	Shoulder Strap	Approx. 376 x 320 x 7.2 mm
	Belt (Sold Separately)	Approx. 1150 x 117 x 8 mm (including extension belt)
	Signal Device	Approx. 42 x 42 x 15 mm
Mass	Main Unit	Approx. 39 g
	Shoulder Strap	Approx. 82 g
	Belt (Sold Separately)	Approx. 51 g (including extension belt)
	Signal Device	Approx. 19 g
Length of Dedicated USB Cable		Approx. 65 cm
Manufacturer		Panasonic Manufacturing (Shanghai) Co., Ltd. No. 258, Jiangtian East Road, Songjiang Industrial Park, Shanghai, China
Country of Origin		China

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© Panasonic Corporation 2026
Panasonic Corporation